

Quick Read 2-in-1 Thermometer

User's Guide TH094

Please be sure to read all warnings, instructions and notations before using this product. Please keep these instructions as they contain important information.

Basics of Your Thermometer

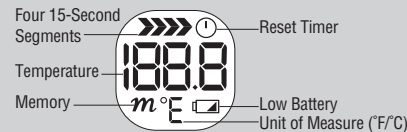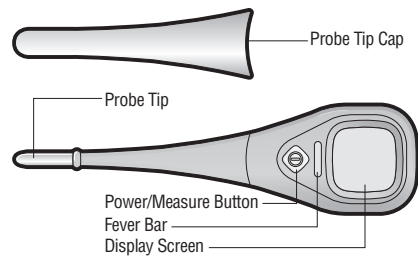Care and Maintenance

Cleaning Probe Tip:

- Clean before and after each use.
- Use an alcohol swab to wipe clean.

Cleaning Thermometer Body:

- Use a soft, dry cloth to clean thermometer body.

Notes:

- **DO** remove battery if stored for a long period of time.
- **DO** store thermometer in a dry place at room temperature with the protective cap in place.
- **DO** take special care not to scratch the display window.
- **DO NOT** allow electronics to get wet.
- **DO NOT** use abrasive cleaning agents, thinners or benzene for cleaning.
- **DO NOT** disassemble for any reason other than battery replacement.
- **DO NOT** drop thermometer.

Normal Body Temperature

Temperature readings are different from person to person and also vary depending upon the body location from which the temperature is measured. For example: Oral temperatures tend to be comparable to ear temperatures for the same individual. Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.

“Normal” temperature refers to each individual’s average temperature while healthy. You may want to measure and record your child’s temperature a few times when he or she is well, so that you will know what is normal for your own child. On average, temperature ranges are as follows:

TEMPERATURE METHOD	NORMAL TEMPERATURE	FEVER RANGE
Oral	97.0° - 99.6°F	99.6°F +
Underarm	95.4° - 98.6°F	98.6°F +

! WARNING BEFORE USING DIGITAL THERMOMETER PLEASE READ CAREFULLY.

• This is **NOT** an ear thermometer. **NEVER** put probe into ear. • This is **NOT** a rectal thermometer. • Take temperature measurements Underarm or Orally. • **ALWAYS** supervise your child when using thermometer. **NEVER** allow your child to use it alone. • **NEVER** allow children to play with this product. It is not a toy. This is a sensitive electronic device and may be damaged if not used properly under parental supervision. • Do not expose the thermometer to extreme temperatures, humidity, direct sunlight. Avoid sudden impact and dust. • Keep battery out of child’s reach. Swallowing a battery could be fatal. • Battery should not be charged or placed into extreme heat, as it may explode. • Used batteries should be disposed of properly.

IMPORTANT INFORMATION • Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique. • **DO NOT** allow child to walk or run during temperature taking. Probe could present a hazard. • High, prolonged fever requires medical attention, especially for young children. Please contact your physician. • Fever chart included on package and back of thermometer.

How to Measure Temperature

Ensure that thermometer is at room temperature, waiting 15 minutes if necessary.

READY DISPLAY:


- Remove probe tip cap. Press and release Power/Measure button. Thermometer beeps once, Fever Bar illuminates amber and all segments display briefly as the thermometer powers on.
- The last temperature measured displays with the *m* icon on screen for a few seconds.
- When Ready Display appears (°F and --- (or °C)) this indicates that the thermometer is ready to take measurement.

Specifics for each Temperature Mode

ORAL—For children 4 years and older.

1. Ensure that °F (or °C) is still flashing, then place probe tip under the tongue in locations indicated by ✓. Instruct children not to bite down on the probe.

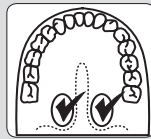- When °F (or °C) stops flashing, and thermometer beeps, the measuring is complete. Approximate measurement time is 8 seconds.
- Remove probe tip from measurement site and check temperature.
- Fever Bar illuminates amber for any temperature above 99.6 °F (37.5 °C), regardless of whether site of measurement is oral or underarm. Consult Fever Chart to note whether the site of measurement indicates an actual fever.

RESET TIMER DISPLAY:

- After reading is complete, the Reset Timer Ⓛ will appear on the screen and begin counting to 60 seconds, showing ▶ every 15 seconds. This allows the thermometer to cool back to room temperature for an accurate second reading.
- Thermometer must complete the Reset Timer Ⓛ cycle to return to room temperature before another measurement can be taken. When Ready Display --- appears this indicates that the thermometer is ready to take measurement.
- To Power off thermometer, press Power/Measure button twice during Reset Timer cycle.

UNDERARM—For child under age 4.

1. Wipe underarm with a dry towel.
2. Ensure that °F (or °C) is still flashing, then place probe tip at center of underarm. Bring arm down against the body. Probe must touch skin.

FEVER CHART	
Underarm	98.6°F + 37°C +
Oral	99.6°F + 37.5°C +

Troubleshooting

Display	Meaning	Recommended Action	Display	Recommended Action
H	<ul style="list-style-type: none"> • Thermometer may have been stored in a hot environment. • Measured temperature is above this thermometer’s range. 	<ul style="list-style-type: none"> • Wait 30 minutes to allow thermometer to reach room temperature, then try again. • See Technical Specifications section for Measurement Temperature Range. 		<ul style="list-style-type: none"> • Batteries are low. Replace them. (See How to Replace the Batteries section).
Err	<ul style="list-style-type: none"> • Thermometer may have been mishandled. You may have irreversibly damaged it. 	<ul style="list-style-type: none"> • Shut off and restart before discarding. 		<ul style="list-style-type: none"> • Do not return this product to the place of purchase. If any parts are missing, email support@djgusa.zendesk.com, or chat with us at www.safety1st.com. You can also call Consumer Care at (800) 544-1108. You can also visit our website at www.safety1st.com. Have ready the model number (TH094) and date code (manufacture date located on package). • Install new battery and re-try function before calling Consumer Relations on electronics.

Memory Recall

Thermometer stores the last temperature reading. As a reference for parents, that temperature reading displays automatically the next time the thermometer is powered on.

Fahrenheit / Celsius

While the thermometer is off, press and hold the Power/Measure button for about 3 seconds. You will see F or a C display on the screen and then hear a short beep. Your thermometer is now in the mode that just appeared on the screen. To switch back, leave thermometer off and hold the button for another 3 seconds.

How to Replace the Batteries


Dead Battery Indicator:

When the “” icon flashes, the battery is exhausted and needs to be replaced.

NOTE: **DO NOT** press or hold down power button at any point while changing the battery or replacing the battery door.

To Replace Battery:

- 1 Use a miniature Phillips head screwdriver to loosen battery cover screw. Remove cover and old battery.
- 2 Insert a new 3 volt battery (CR2032 or equivalent), correctly positioned with the **positive side facing up** as shown in the illustration.
- 3 Place the battery cover on the thermometer and secure the cover.

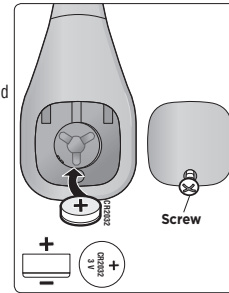Technical Specifications

Measurement Temperature Range: 89.6° F to 109.2° F

Measurement Accuracy: +/- 0.3°F for temp 97.0°F
+/- 0.2°F for temp 101°F
+/- 0.3°F for temp 104°F

Operation Environment: 50° F to 104° F

Storage Environment: -13° F to 140° F

Test Protocol:

- 1 Put the tip of thermometer into water bath, approximately 1/2” below the surface.
- 2 Wait for 15 seconds, then turn on the thermometer by pressing the Power/Measure button. The thermometer is now in standard mode.
- 3 When the measurement is completed, record the temperature reading.

One (1) Year Limited Warranty

Dorel Juvenile Group, Inc. warrants to the original purchaser that this product (Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc., will repair or replace the product, at our option, free of charge. Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on the front of this document and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse.

This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or consequential damages.

SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS THAT MAY VARY FROM STATE TO STATE.

©2017 Dorel Juvenile. (800) 544-1108 www.safety1st.com
Made in CHINA. Hecho en CHINA. Styles and colors may vary. Los estilos y los colores pueden variar. Distributed by (distribuido por) Dorel Juvenile Group, Inc., 2525 State St., Columbus, IN 47201-7494 Dorel Juvenile Canada, 2855 Argentinia Road, Unit 4, Mississauga, ON, L5N 8G6

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