DIGITAL HEATING PAD

ALMOHADILLA CALEFACTORA DIGITAL

PAD CHAUFFANT NUMÉRIQUE

ITEM # RTL19G007GR

INSTRUCTION GUIDE
GUÍA DE INSTRUCCIONES
GUIDE D'INSTRUCTIONS

LEGEND

120 V ~	Alternate current
	Insulation class II
	Instruction for use
	Not to be used by small children (0-3 years)
	Pad is machine washable and can be tumble dried. Must remove plug before washing/drying.
	Do not use pins

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

DANGER- to reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions:

1. Read all instructions.
2. **ALWAYS UNPLUG THE UNIT FROM THE ELECTRICAL OUTLET IMMEDIATELY AFTER USE.**
3. Do not use while sleeping.
4. Do not use on an infant.
5. This pad is not to be used on or by any individuals with certain disabilities, sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
6. Do not use on areas of sensitive skin.
7. Burns can occur regardless of control setting, check skin under pad frequently.
8. Do not use pins or other metallic means to fasten this pad in place.
9. Do not sit on or crush pad - avoid sharp folds.
10. Never pull this pad by the supply cord.
11. Do not use the cord as a handle.
12. Use this pad only on a 110-120 volt AC circuit. Unplug when not in use.
13. Do not tamper with this pad in any way. There are no user serviceable parts.
14. Do not use the pad with liniment, salve or ointment preparations that contain heat-producing ingredients as skin burns could result.
15. Do not use in an oxygen atmosphere.
16. Loop cord loosely when storing, tight wrapping may damage cord and internal parts.
17. Children should be supervised to assure that they do not play with the device.
18. Store in a dry place that is not accessible to children.
19. Do not store until the device returns to ambient temperature.
20. Do not place objects that are either heavy or sharp on top of the device.
21. Use a dry cloth to clean digital remote.
22. The device should not be used by people (including children) that have reduced physical, mental or sensory capabilities nor by people that lack experience or knowledge, unless they are supervised and trained.
23. If the electrical cord is or appears to be damaged, it shall only be replaced by the manufacturer, by the authorized service representative or by a qualified technician in order to eliminate harm or danger.

SAVE THESE INSTRUCTIONS

HEATING PAD COMPONENTS

- Washable Pad with Pad Connector
- Control Connector

CAUTIONS

- This heating pad should not be used in the following cases: if you have diabetes, tuberculosis, benign or malignant tumors, phlebitis and thrombosis, hemorrhages, inflammations of the skin, undetermined calf pain or joint pain, or if you have poor circulation. Consult a physician before using if you have any kind of fungal infection, swelling, fracture or persistent pain.
- **Consult your physician should you have any therapeutic questions regarding use of this heating pad.**
- Do not use on any unexplained pain, or swollen muscles, or following a serious muscular injury before consulting a physician.
- Keep away from children.
- Discontinue use of product if pain increases or swelling occurs after application.
- Do not use while bathing or showering. Use only in a dry area.
- Consult your doctor before using in case of illness or medical condition, including but not limited to pacemaker use, pregnancy or cancer.
- If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.
- Never use any heating products on open wounds, discolored areas, or any area of the body that is swollen, burned, bruised, inflamed or where skin eruptions, ulcerated sores or varicose veins are present.
- This product should never be used by any individual suffering from any ailments that would limit or altogether inhibit the user's capacity to feel or have any sensation in the body part that this product is being used for.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.
- Ensure that your hands are dry when operating the controls or removing the plug.


FCC CAUTION:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

HOW TO ASSEMBLE


The Control Connector should be securely attached to the Pad Connector. Push both parts together for a tight fit.

NOTE: control must be connected to heating pad before plugging unit into outlet.

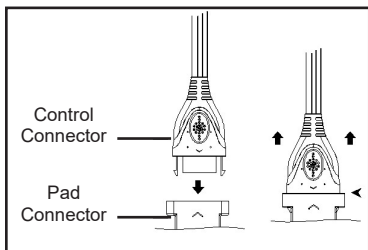TO OPERATE

1. After control is connected to pad, plug into 120V AC electrical outlet. The LCD display will light up and then go out immediately.
2. To turn the pad on, press the Power Button. The LCD display will light up. The default temperature setting is the Lowest setting (113F) and the default timer setting is 1h. This heating pad features fast heat. With this unique fast heat technology, you will feel heat in 30 seconds and can see the real time temperature on the LCD display.
3. Press Temp Button to select temperature settings. The setting grid increases one by one each time you press the button. The temperature figure changes accordingly (113; 122; 131; 140; 149; 158F). The pad starts heating immediately when you release the button. When the pad reaches the desired temperature, it will remain constant for the duration of the time.
4. Press the Timer Button to select Auto-Off timing. The default setting is 1 hour. The 6 timing sequences are: 1:00; 1:30; 2:00; 00:15; 00:30; 00:45. When you release the button, the countdown timing automatically begins.
5. Press Temp Button for 2 to 3 seconds to switch between Fahrenheit and Celsius.
6. The pad and LCD display will shut off automatically when you reach the desired time limit.
7. DURING USE, NEVER PLACE THE CONTROLLER UNDER THE PAD.

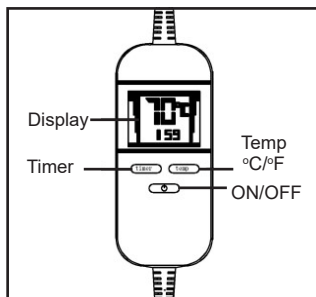MOIST HEAT TREATMENT

Dampen either side of the heating pad with water by patting the surface with a wet sponge or cloth or use a spray bottle to apply a fine mist. Do not thoroughly soak pad.

When using moist heat, you should protect furniture and clothing from becoming damp.

Do not spray or dampen the pad connector.

DANGER: Do not immerse the pad or wiring in water or place under running water. Failure to comply may result in product failure and/or personal injury.

NOTE: The electrical pad connector must be thoroughly dry before establishing electrical connection.

WASHING INSTRUCTIONS

DO NOT DRY CLEAN! DO NOT BLEACH!

DO NOT USE WRINGER FOR DETACHABLE PAD!

1. Disconnect pad from outlet.
2. Detach Pad Connector from the Control Connector.
3. Wash by hand or as follows in an automatic washer using slow agitation and normal spin speed: Fill washer with warm water. Add minimum amount of all-purpose detergent, agitate to dissolve, measure heating pad before washing. Add heating pad and soak five minutes. Agitate two minutes. Spin. Fill for cool rinse, agitate one minute then spin. Dry in automatic dryer on high heat, five minutes only, or line dry. Remove damp heating pad from dryer, block to original size, and drape to finish drying. Do not iron your heating pad. Do not plug in pad and pad connector until they are completely dry.

3 YEAR LIMITED WARRANTY

Your Drive brand product is warranted to be free of defects in materials and workmanship for 3 years of the original consumer purchaser.

This device was built to exacting standards and carefully inspected prior to shipment.

This 3 Year Limited Warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement.

If you have a question about your Drive device or this warranty, please contact an authorized Drive dealer.

© 2021 Medical Depot, Inc. All rights reserved.

Drive is a trademark of Medical Depot, Inc.

Port Washington N.Y. 11050 USA Made in China

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© 2021 Medical Depot, Inc. Tous droits réservés.
Drive est une marque de commerce de Medical Depot, Inc.