**Sick Day Survival Kit**

IH509

Please be sure to read all warnings, instructions and notations before using this product. Please keep these instructions as they contain important information.

**BoogEase (Nasal Aspirator)**

**Parts:**

- (A) Reservoir cap
- (B) Reservoir
- (C) Nose tips
- (D) Mouth piece

**To Assemble:**


**1 Select Nose Tip:** Choose the nose tip that best fits slightly inside your child's nose and creates a seal.

**NOTE:** Aspirator tip is not meant to go deep inside nose, it sits comfortably in nostril opening to create seal.

**2 Tube and Mouth Piece:** Twist and push the tube and mouth piece together and onto to the reservoir of the BoogEase if detached.

**To Use:**  
**1 Position Mouth Piece:** Put the mouth piece in between your lips to prepare for use. Close mouth around mouth-piece and gently hold with front teeth.

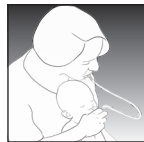**2 Position Child:** Hold your baby facing forward or sit them upright in your lap facing outwards in the crook of one arm.

**3 Position Aspirator:** With your free hand, position the nose piece against their nostril, creating a seal. NOTE: For a more secure grip, it may help to place your index finger between the nose piece and the tube while remaining fingers hold the BoogEase reservoir.

**4 Adult Creates Suction:** Suck air in to create suction in child's nostril. You will be able to see the first chamber begin to fill.

**5 Stopping:** Once the child is decongested or mucus enters the second chamber of the BoogEase, stop inhaling and remove the nose tip from your child's nose.

**To Clean:**

Disassemble as described below to clean parts. Wash with warm, soapy water before first use and after each use. Rinse thoroughly. **DO NOT** microwave.

**Empty:** Separate top and bottom of BoogEase reservoir and empty chambers.

**Aspirator Tube:** Grab the mouth piece at the end of the tube and twist and pull to separate. Repeat action with opposite end of tube. Wet tube brush, then insert into tube to clean inside. Rinse tube in warm soapy water.

NOTE: The cleaning brush is used for both the Aspirator and Fuss Free Medicine Spoon.

**Nose Tip and Reservoir:** Pull nose piece off top of BoogEase reservoir. Wash each piece with warm soapy water.

NOTE: Use provided brush to clean small openings.

**Drying:** Dry with a towel or lay out separated pieces to air dry. Once completely dry reassemble parts.

**⚠ WARNING: CHOKING HAZARD –** Product contains small parts. Keep small parts away from children.

- **NEVER** push tip too far into the nostril and always clean thoroughly after use.
- **NEVER** use in baby's mouth or ears.
- When storing, keep away from heat and bright light.

**Medicine Tracker**

**Purpose:** Hook the top of the medicine tracker around the neck of a medicine bottle and use the bubbles to keep track of when doses are administered. The left column indicates morning and the right columns indicates evening. There are 4 rows, each indicating a new day.

**To Use:**

**Attach:** Hook the top of the medicine tracker around the neck of a medicine bottle. Keep medicine and medicine tracker out of reach of children.

**Mark First Dose:** When you have given your child their first dose, depress bubble correlating with the time of day (morning, or evening).

**Continue Dosage:** For each following dose given that day depress the corresponding bubble.

**New Day:** For each first dose of a day, start on a new row.

**Day 5:** Once you have reached day 5 of doses either:

Pop each bubble back outwards and start on the top row or take tracker off bottle, flip to opposite side, reattach to bottle and start on the top row.

**To Clean:** Wash in warm, soapy water. Product is also top rack dishwasher safe in a basket.

**Fuss Free Medicine Spoon with Syringe and Cleaning Brush**

Makes giving medicine easier by masking its yucky taste with baby's favorite food.

**To Use:**

- 1 Remove syringe from medicine spoon.
- 2 Fill to desired dosage.
- 3 Snap syringe onto medicine spoon.
- 4 Place baby's favorite food on spoon and push plunger slowly to release medicine. Repeat as needed until syringe is empty.

**To Clean:** Wash all parts separately in warm soapy water and rinse thoroughly before first use and after each use. Use provided brush to clean small openings. Product is also top rack dishwasher safe in a basket.

**⚠ WARNING:** Keep all medications out of children's reach. We recommend that you use a child resistant top to seal the medicine bottle and store in a locked cabinet or other location. **MADE WITHOUT BPA.**

**Mesh Bag**

Use mesh bag to keep kit components in. Surface wash only.

**Natural Boogie Wipes**

The moist wipes are provided to help keep baby clean. Lift one end of the cover, take one wipe from the packet and press the cover to seal it, so that the wipes stay moist. Discard wipe after use.

**FOR EXTERNAL USE ONLY.**

**⚠ WARNING:** These products are not toys. **DO NOT** allow children to play with them. To prevent accidental injury, your child should **NEVER** be left unattended. **KEEP OUT OF REACH OF CHILDREN.** We recommend that you store this kit in a locked cabinet or other location that is not accessible to a child. **FOR INDOOR HOME USE ONLY.**

**3-in-1 Digital Thermometer**

The Safety 1st Digital Thermometer takes a clear and accurate temperature without the worries of a glass/mercury device. Effective for oral, rectal or underarm use.

TH093

**Fahrenheit / Celsius**

While thermometer is off, press and hold the ON/OFF button until the unit of measure in the display screen switches. This will take approximately 5 seconds.

**How to Measure Temperature**

Ensure that thermometer is at room temperature, waiting 15 minutes if necessary.

- Press and release ON/OFF button. All segments display briefly as the thermometer powers on.
- The last temperature measured displays with the **M** icon on screen for a few seconds. Then the thermometer will show an average temperature (37.0.9°C or 98.5.9°F) for 0.5 seconds before entering ready mode.
- **Lo** displays and °F flashes when the thermometer is ready to take measurement.

**Specifics for each Temperature Mode**

**ORAL**—For children 4 years and older.

1. Ensure that °F is still flashing, then place probe tip under the tongue in locations indicated by ✓. Instruct children not to bite down on the probe.

**RECTAL**—Highly accurate for measuring core body temperature. Pediatrician recommended mode from birth to age 1, and beyond as tolerated by child.

1. Lubricate probe tip with a water-soluble lubricant or petroleum jelly. For easier insertion, position baby so that the back and legs form a 90-degree angle (L-shape).
2. Ensure that °F is still flashing, then insert thermometer tip gently into the rectum, taking care not to insert past the insertion gauge (ring on the probe tip). Do not force the tip into rectum if resistance is encountered.

- When °F stops flashing, the measuring is complete. Approximate measurement times are: 60 seconds for oral and rectal temps, 90 seconds or longer for underarm temp.
- Remove probe tip from measurement site. Read and record the temperature and time for reference.
- Press button to shut off or thermometer will shut off automatically after 8-9 minutes.

NOTE: High, prolonged fever requires medical attention, especially for young children. Please contact your physician. Fever chart included. **DO NOT** allow child to walk or run during temperature taking. Probe could present a hazard.

**Normal Body Temperature**

Temperature readings are different from person to person and also vary depending upon the body location from which the temperature is measured. For example: Oral temperatures tend to be comparable to ear temperatures and approximately 1° lower than rectal temperatures for the same individual. Temperature readings are affected by many factors, including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.

"Normal" temperature refers to each individual's average temperature while healthy. You may want to measure and record your child's temperature a few times when he or she is well, so that you will know what is normal for your own child. On average, temperature ranges are as follows:

TEMPERATURE METHOD	NORMAL TEMPERATURE	FEVER RANGE
Oral	97.0° - 99.6°F	99.6°F +
Rectal	97.3° - 100.4°F	100.4°F +
Underarm	95.4° - 98.6°F	98.6°F +

**One (1) Year Limited Warranty**

Dorel Juvenile Group, Inc. warrants to the original purchaser that this product (Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc., will repair or replace the product, at our option, free of charge. Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on the front of this document and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse.

This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or consequential damages.

**SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.**

**THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS THAT MAY VARY FROM STATE TO STATE.**

**⚠ WARNING:** Before using digital thermometer, please read carefully.

- This is **NOT** an ear thermometer. **NEVER** put probe into ear.
- Take temperature measurements Underarm, Orally or Rectally.
- **ALWAYS** supervise your child when using thermometer and thermometer covers.
- **NEVER** allow your child to use it alone or without close parental supervision.
- **NEVER** allow children to play with this product. It is not a toy. This is a sensitive electronic device and may be damaged if not used properly under parental supervision.
- **DO NOT** expose the thermometer to extreme temperatures, humidity, direct sunlight or shock. Avoid sudden impact and dust.
- Keep battery out of child's reach. Swallowing a battery could be fatal. Battery should not be charged or placed into extreme heat as it may explode.
- Used batteries should be disposed of properly.

**Technical Specifications**

**Measurement Range:** 90° F to 109.9° F

**Measurement Accuracy:** +/- 0.2° F for temps 96° F to 107.6° F  
+/- 0.4° F for temps below 96° F or above 107.6° F.

**Operation Environment:** 60.8° F to 104° F

**Storage Environment:** -13° F to 144° F

**Test Protocol:**

- 1 Put the tip of thermometer into water bath, approximately 2 cm below the surface.
- 2 Wait for 10 seconds, then turn on the thermometer by pressing the ON/OFF button. The thermometer is now in standard mode.
- 3 When the measurement is completed (approximately 60 seconds) record the temperature reading.

**UNDERARM**—For child under age 4 who will not cooperate with rectal temperature measurement.

1. Wipe underarm with a dry towel.
2. Ensure that °F is still flashing, then place probe tip at center of underarm. Bring arm down against the body. Probe must touch skin.

**How to Replace the Batteries**

**Dead Battery Indicator:**

- When the "□" icon appears, the battery is exhausted and needs to be replaced.
- Battery is to be inserted with the correct polarity.
- Check that all contact surfaces are clean and bright before installing battery.

**To Replace Battery:**

- 1 Use a Phillips head screwdriver to loosen battery cover screw. Remove cover.
- 2 Use a toothpick to push the battery out of the metal housing.
- 3 Insert a new battery (LR41 or equivalent), correctly positioned with the positive "+" side facing up as shown in illustration.
- 4 Place the battery cover on the thermometer and secure the cover.

**Care and Maintenance**

**Cleaning Probe Tip:**

- Clean before and after each use.
- Use an alcohol swab to wipe clean.

**Cleaning Thermometer Body:**

- Use a soft, dry cloth to clean thermometer body.

**Notes:**

- **DO** remove battery if stored for a long period of time.
- **DO** store thermometer in a dry place at room temperature with the protective cap in place.
- **DO** take special care not to scratch the display window.
- **DO NOT** allow electronics to get wet.
- **DO NOT** use abrasive cleaning agents, thinners or benzene for cleaning.
- **DO NOT** drop thermometer.
- Install new battery and retry function before calling Consumer Relations on electronics.

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*Please be sure to read all warnings, instructions and notations before using this product. Please keep these instructions as they contain important information.*

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