Gentle Read Rectal Thermometer

TH013 & TH035 User Guide

! WARNING: BEFORE USING DIGITAL THERMOMETER, PLEASE READ CAREFULLY.

- This is **NOT** an ear thermometer. **NEVER** put probe into an ear.
- This is a rectal thermometer **ONLY**. It should **NEVER** be used for oral, underarm, ear or any other temperature measurement.
- **NEVER** allow children to play with this product. It is not a toy. This is a sensitive electronic device and may be damaged if not used properly under parental supervision.
- **ALWAYS** supervise your child when using thermometer. **NEVER** allow your child to use it alone or without close parental supervision.
- Avoid sudden impact and dust.
- **DO NOT** expose the thermometer to extreme temperatures, humidity, direct sunlight or shock.


Know Your Unit

Product Features:

- Fast, 8-second rectal readings.
- Memory recall of last reading.
- Long life, replaceable battery (approximately 2000 temperatures).

For Product Questions or Problems:

- Visit www.safety1st.com and look for "Safety 1st TV" for information about how to take a rectal temperature.
- Do not return this product to the place of purchase. If any parts are missing, email consumer@djgusa.com, call Consumer Relations at (800) 544-1108, or fax at (800) 207-8182. You can also visit our website at www.safety1st.com and follow links to "Contact Us." Have ready the model number (TH013 or TH035) and date code (manufacture date located on package and in battery compartment).
- Install new battery and retry function before calling Consumer Relations on electronics.

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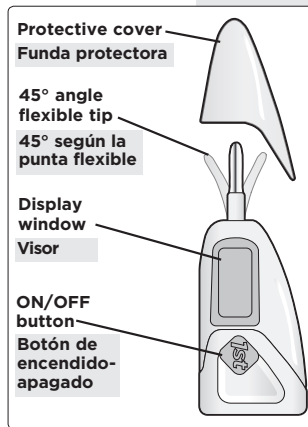
To Change Between Fahrenheit and Celsius

Ensure that thermometer is off. Then press and hold the ON/OFF button for at least 2 seconds. When you stop pressing the ON/OFF button after 2 seconds, the current measurement scale (°F or °C icon) will display. Change the measurement scale between °F and °C by pressing the ON/OFF button again. When the measurement scale has been chosen, wait for 3 seconds, and the unit will enter the "Ready for Measuring" mode automatically.

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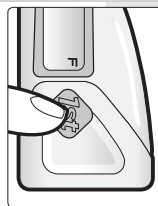
To Replace Button Cell Battery (LR41)

When the "▽" appears in the lower right of the display, the battery is exhausted and needs replacing. Make sure a battery of the same type is on hand and proceed as follows after locating the battery compartment cover on the back side of the thermometer.

- 1 Carefully unscrew battery compartment cover using a Phillips head screwdriver. Slide cover downward to remove it.
- 2 Use a toothpick to remove battery. Insert new equivalent battery type (LR41). Ensure battery is correctly positioned with the "-" facing upward as shown in the illustration. Carefully replace the battery compartment cover and secure the screw. Properly dispose of the batteries, keeping them from small children and heat.

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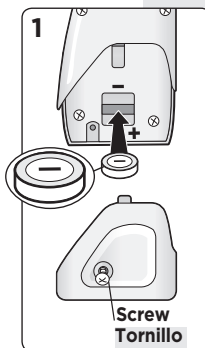! WARNING:

- Keep battery out of child's reach. Swallowing battery could be fatal.
- Battery should not be charged or placed into extreme heat as it may rupture.
- Used batteries should be disposed of properly.

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How to Take Temperature

Commonly used for babies and young children or when it is difficult to take an oral or underarm temperature.

NOTE: Be sure probe tip is clean before using. See Care and Cleaning.

- 1 Push the ON/OFF button located beside the display window. The display will read:

This is a function check indicating the unit is working properly. The last measured temperature will be displayed on the screen for two seconds. The thermometer will read "L" and begin measuring temperature when exposed to body heat.

- 2 The °F symbol will continuously flash until after completion of the temperature taking process. Please note that the "▼" indicates a low battery. See instructions for replacing battery before measuring temperature. An "ERR" message in the window indicates your unit is not functioning properly. Do not expect an accurate temperature reading if this low battery indicator is present. Out-of-range indicator displays "L" when measurement is lower than 90.0°F and "H" when higher than 109.0°F.
- 3 Before using, lubricate probe tip with a water-soluble lubricant or petroleum jelly. For easier insertion, position baby so his back and legs form a 90-degree (L-shaped) angle. Try placing child on your lap with his legs down, or lay him on a bed or changing table using a small pillow to raise his hips slightly. Spread his buttocks with one hand while using the other to gently insert the thermometer tip into the rectum, **never more than 1/2 inch**. The short probe tip ensures you don't insert too far. **Never force the thermometer**. Hold it gently in place until reading is complete. The peak temperature should be reached in approximately 8 seconds.
- 4 The thermometer will beep repeatedly and the °F symbol will stop flashing following peak temperature reading. **Please note: Peak temperature is reached in 8 seconds. Audible beeps may not always be heard. Remove thermometer after 8 seconds. If the thermometer stays in the measuring site, the temperature may rise slightly over time even when the °F does not flash.**
- 5 Read and record the temperature and time for reference. The last reading will automatically be stored in the memory. **For a rectal temperature, a reading of 100.4 °F or higher indicates a fever.**
- 6 Push the ON/OFF button to turn the unit OFF. If you do not turn the unit OFF, it will SHUT OFF automatically in approximately 10 minutes.

Technical Specifications

Measurement Range: 90.0°F to 109.0°F.

Measurement Accuracy: 0.2°F between 96.0°F to 107.0°F.

Test Protocol:

- 1 Put the tip of thermometer into water bath, approximately 1/2" below the surface.
- 2 Wait for 10 seconds, then turn on the thermometer by pressing the ON/OFF button. The thermometer is now in standard mode.
- 3 When the measurement is completed, record the temperature reading.

Care and Cleaning

Do:

- Clean the thermometer before first use and after every use
- Clean thermometer probe tip (1 inch only) using soap and water, or sterilize probe tip (1 inch only) using isopropyl alcohol or other sterilizing solution
- Clean plastic parts of this thermometer with a soft, dry cloth
- Store thermometer in its protective cover when not in use
- Examine thermometer before each use for cracks or other signs of wear
- Discard a worn or damaged thermometer immediately

Do Not:

- Never immerse electronic parts in water, isopropyl alcohol or any other liquid
- Never attempt to sterilize thermometer in boiling water
- Do not store or keep in direct sunlight or in high temperatures
- Do not bend probe tip beyond 45° (to avoid damage to the probe)
- Do not attempt to pull or twist the probe apart from the body or to detach the "1st" button head

One (1) Year Limited Warranty

Dorel Juvenile Group, Inc., warrants to the original purchaser that this product (Gentle Read Rectal Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc., will repair or replace the product, at our option, free of charge. Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on the front of this document and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse.

This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or consequential damages.

SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS THAT MAY VARY FROM STATE TO STATE.

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