
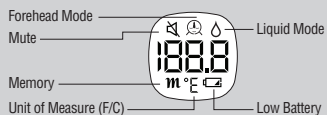
## Easy Read Dual Mode Forehead Thermometer

User's Guide TH091

Please be sure to read all warnings, instructions and notations before using this product. Please keep these instructions as they contain important information.

### Basics of Your Thermometer
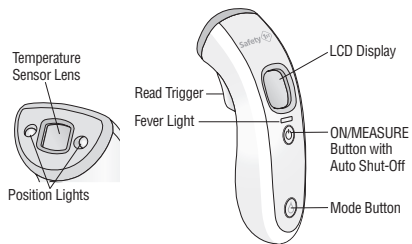
### Care and Maintenance

#### Cleaning Temperature Sensor Lens:

- Clean very gently with a soft, alcohol-moistened cloth.
- Do not use an excessive amount of alcohol.
- Do not pour alcohol directly on sensor.

#### Cleaning Thermometer Body:

- Use a soft, dry cloth to clean thermometer body.

#### Notes:

- **DO** remove battery if stored for a long period of time.
- **DO** store thermometer in a dry place at room temperature.
- **DO** take special care not to scratch the sensor or display window.
- **DO NOT** use abrasive cleaning agents, thinners or benzene for cleaning.
- **DO NOT** expose thermometer to extreme temperatures, humidity or direct sunlight.
- **DO NOT** drop thermometer.
- **DO NOT** submerge any part of the product in water.

**PLEASE NOTE:**  
The sensor lens is the most delicate part of the thermometer. To ensure accurate readings it must be clean and intact.

### Normal Body Temperature

Temperature readings are different from person to person and also vary depending upon the body location from which the temperature is measured. Forehead temperatures tend to be comparable to oral temperatures and approximately 1° lower than rectal temperatures for the same individual. Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.

"Normal" temperature refers to each individual's average temperature while healthy. You may want to measure and record your child's temperature a few times when he or she is well, so that you will know what is normal for your own child. On average, temperature ranges are as follows:

TEMPERATURE METHOD	NORMAL TEMPERATURE	FEVER RANGE
Forehead	97°F to 99.5°F	99.6°F +

normal forehead temp range

### Mute Function

To set to mute: While unit is ON, press the ON/MEM button for 3 seconds until you see the "M" icon appear. Then release ON/MEM button.

### Memory

While unit is ON, press the ON/MEM button. The number "1" will appear with the "M" icon to indicate this is the first reading in the memory. Release the ON/MEM button to see the reading. Press the ON/MEM button again to toggle through up to 25 records.

**WARNING** • This is **NOT** an ear thermometer. Take temperature measurements from the center of forehead. • Use caution when taking temperature readings of bath water and food. This product is not intended to touch the person or item whose temperature is being measured. • **ALWAYS** supervise your child when using the thermometer. **NEVER** allow your child to use it alone. • **NEVER** allow children to play with this thermometer. It is not a toy. This is a sensitive electronic device and may be damaged if not used properly under parental supervision. • **DO NOT** expose the thermometer to extreme temperatures, humidity, direct sunlight or shock, and/or dust. • Keep battery out of child's reach. Swallowing a battery could be fatal. Battery should not be charged or placed into extreme heat as they may explode. • Used batteries should be disposed of properly. **Other Battery Information:** We recommend alkaline batteries for longer battery life. **DO NOT** mix old and new batteries. **DO NOT** mix alkaline, standard (carbon zinc) or rechargeable (nickel-cadmium) batteries. Non-rechargeable batteries are not to be recharged. Rechargeable batteries are to be removed from the unit before being charged. Rechargeable batteries are only to be charged under adult supervision. Only batteries of the same or equivalent type as recommended are to be used. Batteries are to be inserted with the correct polarity. The supply terminals are not to be short-circuited. Check that all contact surfaces are clean and bright before installing batteries. **DO NOT** disassemble for any reason other than inserting or replacing batteries.

### How to Measure Temperature *Make sure the thermometer has been at room temperature for 30 minutes before use.*

#### 1. Turn Thermometer On and Prep Forehead

Press and release the ON/MEMORY button. Thermometer runs a brief self-check, then beeps and displays the ready screen. The unit will default to FOREHEAD MODE. POSITION LIGHTS at front of unit come on automatically. Prep the forehead by making sure it is clear of hair and sweat. If child has recently been in a very warm or very cold environment, wait a few minutes for their skin to adjust to room temperature before taking their temp.

#### 2. Position Thermometer

The POSITION LIGHTS are LED lights and are meant to help you feel confident about how you are aiming and distancing the unit. **-AIMING:** The SENSOR LENS points at the area of the forehead between the two red POSITION LIGHTS. The right spot to aim the SENSOR LENS is at the base of the forehead, just above the eyebrows. So, when properly positioned, the LOWER POSITION LIGHT should be between the eyebrows. **-DISTANCE:** For accurate temperatures, the unit should be within 1.5 inches (4 cm) of the forehead. When the unit is 1.5 inches from the forehead, the UPPER POSITIONING LIGHT will be directly above the LOWER POSITIONING LIGHT (vertically.) When closer or farther away, the lights move away from each other in opposite directions. **Hint:** If a child is uncomfortable seeing the thermometer coming toward his face, try turning the thermometer upside-down. The POSITIONING LIGHTS work the same when the thermometer is upside down.

#### 3. Measure Temperature

Once unit is positioned, press and hold the READ TRIGGER button on the rear of the thermometer for about 1 second or until the screen lights up, displays temperature reading, and beeps (if not in mute mode.)

#### 4. Fever Reading Indicator

The unit gives one long beep when temp is in normal range. Thermometer gives additional beeps and the fever bar lights up when elevated temp is detected. **NOTE: High, prolonged fever requires medical attention, especially for young children. Please contact your physician. Fever chart included.**

#### 5. Wait timer & Auto Shut-Off

After a reading is complete, the WAIT icon will flash for 5 seconds. After those 5 seconds are complete, the unit will return to READY mode. Thermometer beeps and shuts off automatically after about 30 seconds in READY mode.

#### 6. Liquid Mode

This thermometer can also be used to take no-touch temperatures of liquids like bath water and milk/formula. Unit defaults into FOREHEAD MODE every time you turn it on. To switch to LIQUID MODE, simply press the MODE button on the handle of the thermometer. **NOTE:** To return to FOREHEAD MODE, press MODE button again. **NOTE:** Once in LIQUID MODE, positioning lights and fever bar are no longer needed.

**NOTE:** Room temperature is between 50°F and 104°F (10°C - 40°C).

Display	Meaning	Recommended Action
Er	• Error 3, the ambient temperature is not within range • Error 5-9, the system is not functioning properly	• Allow the thermometer to rest in a room for at least 15 minutes at room temperature. • Remove batteries, wait 1 minute and repower unit. If the message reappears, contact Consumer Care.
H i	• Forehead Mode: Temperature taken is higher than 108°F (+42.2°C) • Liquid Mode: Temperature taken is higher than 176°F (+80°C)	• Allow the thermometer to rest in a room for at least 15 minutes at room temperature.
-- --	• Forehead Mode: Temperature taken is lower than 89.6°F (+32°C) • Liquid Mode: Temperature taken is lower than -7.6°F (-22°C)	• Target is too cold to read using this thermometer. If in forehead mode, try moving unit closer to the forehead.
Battery Icon	• Device cannot be powered on to the ready stage	• Remove batteries, wait 1 minute, and repower unit. If the message reappears, contact Consumer Care.

### Troubleshooting

Display	Recommended Action
Battery Icon	• Batteries are low. Replace them. (See Battery Information section).
Wait Icon	• For best accuracy, wait one minute between temperature measurements so that thermometer has time to re-acclimate to the room's temperature.
Memory Icon	• Do not return this product to the place of purchase. If any parts are missing, email support@digiusa.zendesk.com, or chat with us at www.safety1st.com. You can also call Consumer Care at (800) 544-1108. Have ready the model number (TH091) and date code (manufacture date located on package).
Ready Icon	• Install new batteries and re-try function before calling Consumer Care on electronics.

### How to Install/Replace the Batteries

This thermometer comes with the batteries outside of the unit. If this is your first time using this unit, start with step 1 below.


#### Low Battery Warning:

- When the icon appears on the display screen when powering up and using the thermometer, the batteries are low but may still work for a very short time.
- When the icon appears and flashes by itself on the screen the batteries are exhausted and need to be replaced.


#### To Install Batteries:

- 1 Use a Phillips head screwdriver to loosen battery cover screw. Remove cover.
- 2 Remove old batteries
- 3 Insert two AAA (or equivalent) batteries, correctly positioned with positive "+" and negative "-" as shown in illustration.
- 4 Place the battery cover on the thermometer and tighten the screw to secure it in place.
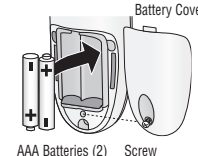
### Switch Between Fahrenheit & Celsius

- 1 If unit is not off already, hold the POWER/MEMORY button for 6 seconds. Release when you see the word "OFF" appear on the screen.
- 2 While unit is still off, hold down READ TRIGGER button and then POWER/MEMORY button. Hold both down for about 5 seconds.
- 3 After about 5 seconds, unit will power on in the new mode.

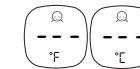**NOTE:** Switching between F/C clears the memory on the unit.

### Technical Specifications

**Measurement Temperature Range:**  
**Forehead Mode:** 89.6°F to 108.0°F (32°C to 42.2°C)  
**Liquid Mode:** -7.6°F to 176.0°F (-22°C to 80°C)

**Measurement Accuracy:**  
**Forehead Mode:** +/- 0.4°F (0.2°C) within 95°F to 107.6°F (35 to 42°C)  
 +/- 0.5°F (0.3°C) for other range  
**Liquid Mode:** +/- 0.5°F (0.3°C) within 71.6°F to 108°F (22 to 42.2°C)  
 others +/- 0.4°F or +/- 0.4°F (0.2°C) whichever is greater.

**Operation Environment:** 50°F to 104.0°F (10 to 40°C)  
 Humidity up to 85% (non condensing)

**Storage Environment:** -4°F to 122°F (-20 to +50°C)  
 Humidity up to 85% (non condensing)

### One (1) Year Limited Warranty

Dorel Juvenile Group, Inc. warrants to the original purchaser that this product (Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc. will repair or replace the product, at our option, free of charge. Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on the front of this document and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse.

This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or consequential damages.

**SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES. SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.**

**THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS THAT MAY VARY FROM STATE TO STATE.**

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©2019 Dorel Juvenile Group, Inc. Safety 1st™ and the Safety 1st logo and design are trademarks of Dorel Juvenile Group, Inc. (800) 544-1108 www.safety1st.com. Styles and colors may vary. Made in CHINA. Distributed by Dorel Juvenile Group, Inc., 2525 State St., Columbus, IN 47201-7494. Dorel Juvenile Canada, 2855 Argente Road, Unit 4, Mississauga, ON L5N 8G6. Importado y Distribuido en Chile por: Comercial e Industrial Sifa S.A. - San Ignacio 0201 / Parque Industrial Portezuelo Quilicura, Santiago Chile. Fono (562) 2339 9000. En Perú por: Comesa S.A. - Los Libertadores 455 - San Isidro, Lima - Perú. Fono (511) 4469014. En Centro América, Caribe y Ecuador por: Best Brands Group S.A. - Av. Balboa, Bay Mall Plaza, Planta Baja, local 9, Panamá - Fono: (507) 300 2884. En Colombia por: Baby Universe S.A.S. - Vía 40 No. 77-29, Barranquilla, Colombia - Fono: (575) 353 1110. Importador: DJGM, S.A. DE C.V. Gabriel Manzano No. 1041 Col. Del Valle, Benito Juárez C.P. 03100 México, Ciudad de México R.F.C. DJG140305GN9 Tel. (55) 67199202 0425/19 4358-7770