Homedics

Shiatsu Air Pro Foot Massager


Register your product today at
www.homedics.co.uk/product-registration

3 YEAR GUARANTEE

FMS-350H-EU

PRODUCT FEATURES

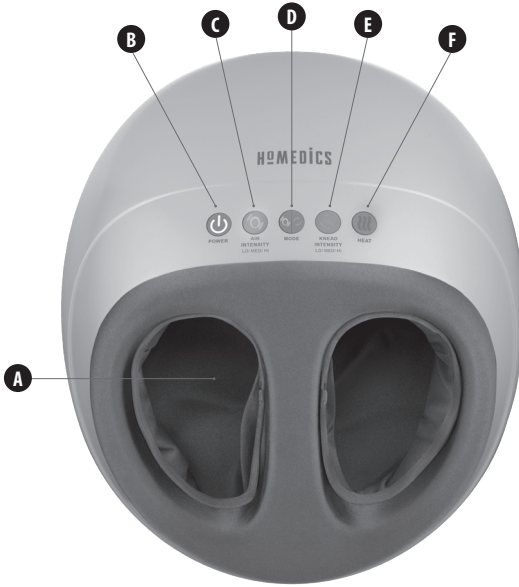- A / ROTATING MESSAGE NODES** - For tired soles and arches
- B / POWER BUTTON** - For Activates massage and heat functions
- C / AIR INTENSITY BUTTON** - Choose low, medium, or high compression intensity
- D / MODE BUTTON** - Activates air only, knead only, or a combination of both
- E / KNEAD INTENSITY BUTTON** - Choose low, medium, or high kneading intensity
- F / HEAT BUTTON** - Soothes tired soles and arches

INSTRUCTIONS FOR USE:

1. Plug into a 220-240V outlet.
2. It is highly recommended that you wear socks while using this foot massager. Do not use the unit with shoes.
3. Place massager on the floor in front of your chair or seat and place feet in massager.
4. To activate massage, press the **POWER button** once. The button will illuminate white to confirm your selection. The unit will start with both air and knead at low intensity.
5. To switch the unit off, press the **POWER button** once.

Air Intensity Button

1. Pressing the **AIR INTENSITY** button will allow you to adjust the intensity of the air pressure, how firmly it squeezes your feet.
2. Press **MODE** button until **AIR INTENSITY** button is illuminated.
3. When first turned on, **AIR INTENSITY** will begin in the lowest setting with the button illuminated dull blue.
4. Pressing the button again will change the foot pressure to medium, indicated by a brighter blue.
5. Pressing the button a third time will change the pressure to high, indicated by the brightest blue.
6. Pressing the button once more will return the pressure to its lowest setting.

Mode Button

1. Pressing the **MODE** button starts the unit with both air and knead in low intensity. Both sides of the button will be illuminated.
2. Pressing the button again will change the unit to air only. The left side of the **MODE** button will illuminate blue.
3. Pressing the button a third time will change the unit to knead only. The right side of the **MODE** button will illuminate green.
4. Pressing the button once more will return the unit to air and knead. **AIR INTENSITY**, **MODE**, and **KNEAD INTENSITY** buttons will all be illuminated.

Knead Intensity Button

1. Pressing the **KNEAD INTENSITY** button will allow you to adjust the intensity of the massaging mechanism.
2. Press **MODE** button until the **KNEAD INTENSITY** button is illuminated.
3. When first turned on, **KNEAD INTENSITY** will begin in the lowest setting with button illuminated dull green.
4. Pressing the button again will change the massage intensity to medium, indicated by a brighter green.
5. Pressing the button a third time will change the intensity to high, indicated by the brightest green.
6. Pressing the button once more will return the intensity to its lowest setting.

Heat Button

1. To activate soothing heat, first select **AIR** or **KNEAD**. Then simply press the **HEAT** button and it will illuminate red.
 2. To turn heat off, press the button again and the red light will turn off.
- NOTE:** For your safety, heat cannot be turned on unless air and/or knead program is selected.

NOTE:

- DO NOT attempt to stand on this Foot Massager. Use only while seated.
- DO NOT use the massage or heat functions continuously for more than 15 minutes.
- There is a 15 minute auto shut-off on this appliance for your safety. This feature should not be considered as a substitute for 'OFF'. Always remember to switch the appliance off when not in use.

MAINTENANCE:

Cleaning

- Unplug the appliance and allow to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Never allow water or any other liquids to come into contact with the appliance.
- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, glass / furniture polish, paint thinner etc to clean.

READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- NEVER leave an appliance unattended when plugged in. Unplug from the mains outlet when not in use and before putting on or taking off parts or attachments.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins or other metallic fasteners into the appliance.

Storage

Place the appliance in its box or in a safe, dry, cool place. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.

- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics.
- NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair.
- Keep cord away from heated surfaces.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the cord or use cord as a handle.
- DO NOT use outdoors.
- DO NOT crush. Avoid sharp folds.
- This product requires a 100-240V AC power.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. All servicing of this appliance must be performed by authorised HoMedics service personnel only.
- Never block the air openings of the appliance. Keep air openings free of lint, hair etc.
- Burns can result from improper use.
- Please ensure that all hair, clothing and jewellery are kept clear of the massage mechanism or any other moving parts of the product at all times.

- If you have any concerns regarding your health, consult a doctor before using this appliance.
- This is a non-professional appliance, designed for personal use and intended to provide a soothing massage to worn muscles. **DO NOT** use as a substitute for medical attention.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your GP.
- Pregnant women, diabetics and individuals with pacemakers should consult a doctor before using this appliance. Not recommended for use by individuals with sensory deficiencies including diabetic neuropathy.
- **DO NOT** use on an infant, invalid or on a sleeping or unconscious person. **DO NOT** use on insensitive skin or on a person with poor blood circulation.
- This appliance should **NEVER** be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- **NEVER** use directly on swollen or inflamed areas or skin eruptions.
- **DO NOT** use before going to bed. The massage has a stimulating effect and can delay sleep.
- Do not use for longer than the recommended time.
- Only gentle force should be exerted against the massage mechanism in order to eliminate risk of injury.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- Failure to follow the above may result in the risk of fire or injury.

3 YEAR GUARANTEE

FKA Brands Ltd guarantees this product from defect in material and workmanship for a period of 3 years from the date of purchase, except as noted below. This FKA Brands Ltd product guarantee does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorised accessory; alteration to the product; or any other conditions whatsoever that are beyond the control of FKA Brands Ltd. This guarantee is effective only if the product is purchased and operated in the UK / EU. A product that requires modification or adaptation to enable it to operate in any country other than the country for which it was designed, manufactured, approved and / or authorised, or repair of products damaged by these modifications is not covered under this guarantee. FKA Brands Ltd shall not be responsible for any type of incidental, consequential or special damages.

To obtain guarantee service on your product, return the product post-paid to your local service centre along with your dated sales receipt (as proof of purchase). Upon receipt, FKA Brands Ltd will repair or replace, as appropriate, your product and return it to you, post-paid. Guarantee is solely through HoMedics Service Centre. Service of this product by anyone other than the HoMedics Service Centre voids the guarantee. This guarantee does not affect your statutory rights.


For your local HoMedics Service Centre, go to
www.homedics.co.uk/servicecentres

WEEE explanation

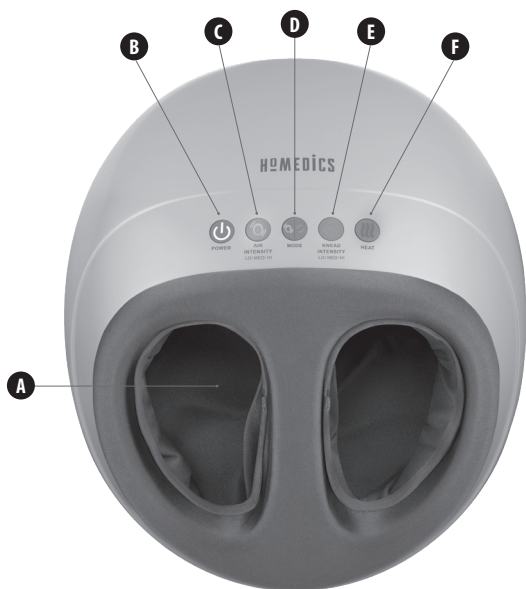This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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FKA Brands Ltd, Somerhill Business Park, Tonbridge, Kent TN11 0GP, UK
FKA Brands Ltd, 29 Earlsfort Terrace, Dublin D02AY28, Ireland