

Homedics

Shiatsu neck massager with heat



✓ Register your product today at
www.homedics.co.uk/product-registration

3 YEAR GUARANTEE

NMS-65H-GB

PRODUCT FEATURES:



Fig.1

Fig.2



The actual product materials and fabrics may vary from the picture shown.

INSTRUCTIONS FOR USE:


1. Plug adaptor into a 220-240V AC outlet and plug the receptacle end into the input jack, located along the handle of the neck massager.
2. Place massager around neck (Fig. 1) or area of back you wish to massage. Convenient straps help assist in easy positioning and intensity control.

Note: Do not use this product for more than 15 mins at a time.

There is a 15 minute auto shut-off on this appliance for your safety. This feature should not be considered as a substitute for 'OFF'. Always remember to switch the appliance off when not in use.

Using essential oils with the aroma tab (Fig. 2)

Gently pour a few drops of oil onto the aroma tab, releasing the scent in the air. The more drops, the stronger the scent.

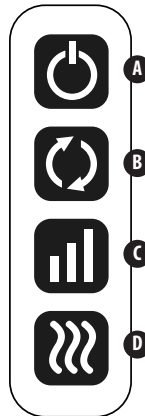

A / POWER BUTTON

To turn on the massage function, press the power button once. Press a second time to switch off.

B / DIRECTION BUTTON

Press the direction button to switch between clockwise and anti-clockwise massage direction.

C / SPEED CONTROL

The massager will start up at its highest speed setting. Press the speed control button during use to change the speed between the 3 settings - high, med and low.

D / HEAT BUTTON

For soothing heat, when massage is on, simply press the heat button and the massage nodes on the neck unit will illuminate. To deselect, press the button again and the corresponding massage nodes will stop glowing.

MAINTENANCE:

Cleaning

Unplug the appliance and allow to cool before cleaning. Clean only with a soft, slightly damp sponge.

- Never allow water or any other liquids to come into contact with the appliance.
- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, glass / furniture polish, paint thinner etc to clean.

Storage

Place the appliance in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.

**READ ALL INSTRUCTIONS BEFORE USE.
SAVE THESE INSTRUCTIONS FOR FUTURE
REFERENCE.**

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- NEVER leave an appliance unattended when plugged in. Unplug from the mains outlet when not in use and before putting on or taking off parts or attachments.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins, metallic fasteners or objects into the appliance or any opening.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics. NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair.
- Keep cord away from heated surfaces.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the cord or use cord as a handle.
- DO NOT use outdoors.
- DO NOT crush. Avoid sharp folds.
- This product requires a 220-240 V AC power.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. All servicing of this appliance must be performed at an authorised HoMedics Service Centre.
- DO NOT sit or stand on the massage part of the appliance. The unit should only be used when placed on the floor and user is seated.
- Never block the air openings of the appliance. Keep air openings free of lint, hair etc.
- Burns can result from improper use.
- Please ensure that all hair, clothing and jewellery are kept clear of the massage mechanism or any other moving parts of the product at all times.
- If you have any concerns regarding your health, consult a doctor before using this appliance.
- This is a non-professional appliance, designed for personal use and intended to provide a soothing massage to worn muscles. DO NOT use as a substitute for medical attention.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your GP.
- Pregnant women, diabetics and individuals

- with pacemakers should consult a doctor before using this appliance. Not recommended for use by individuals with sensory deficiencies including diabetic neuropathy.
- DO NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation.
 - This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
 - NEVER use directly on swollen or inflamed areas or skin eruptions.
 - DO NOT use before going to bed. The massage has a stimulating effect and can delay sleep.
- Do not use for longer than the recommended time.
 - Only gentle force should be exerted against the massage mechanism in order to eliminate risk of injury.
 - The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
 - Failure to follow the above may result in the risk of fire or injury.

Homedics

Distributed in the UK by

FKA Brands Ltd, Somerhill Business Park, Tonbridge, Kent TN11 0GP, UK

EU Importer

FKA Brands Ltd, 29 Earlsfort Terrace, Dublin 2, Ireland

Customer Support: +44(0) 1732 378557 | support@homedics.co.uk

IB-NMS65HGB-0622-01

