Homedics

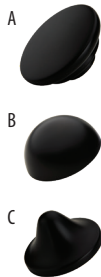
**portable full-body
vibration massager**

Register your product today at
www.homedics.co.uk/product-registration

3 YEAR GUARANTEE

HHP-120BK

INSTRUCTIONS FOR USE

1. This unit comes with two AAA batteries. Put the two AAA batteries into the battery compartment of the unit.
2. To turn the unit on, slide the power  button upward once. This will activate the low intensity vibration setting. To switch to the high intensity vibration setting, slide the power  button up one more time (Fig. 1).
3. Apply the massager head lightly to the area you wish to massage. Move the unit slowly, holding it over the targeted area. Massage the area for a short period of time and move to another area. NEVER massage the same area for longer than 3 minutes.
4. To turn the unit off, slide the power  button all the way down until the vibration stops.

ATTACHMENTS

- A. Flat attachment - moulding muscles all over body
- B. Round attachment - pecs, glutes, quads, lats, hamstrings. suitable for full-body use
- C. Trigger point attachment - suitable for impacting deep tissue and trigger points, neck, palms, soles

MAINTENANCE & STORAGE

Turn the device off and let it cool down before cleaning it. Clean only with a damp cloth and then dry completely. Do not immerse in any liquid such as water or alcohol, steel wool, abrasive cleaners or corrosive liquids to clean it. Place the appliance in its box or in a safe, dry, cool place.

READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

This appliance can be used by children aged from 16 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision. Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid. DO NOT reach for an appliance that has fallen in water or other liquids. Keep dry – DO NOT operate in wet or moist conditions. Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics. NEVER operate the appliance if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair. DO NOT attempt to repair the appliance. There are no user serviceable parts. All servicing of this appliance must be performed at an authorised HoMedics Service Centre. Please ensure that all hair, clothing and jewellery are kept clear of moving parts of the product at all times. If you have any concerns regarding your health, consult a doctor before using this appliance. Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your GP. Pregnant women, diabetics and individuals with pacemakers should consult a doctor before using this appliance. Not recommended for use by individuals with sensory deficiencies including diabetic neuropathy. DO

NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation. This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls. Do not use this product for more than 15 mins at a time. Use this product only on the body's soft tissue as desired without producing pain or discomfort. Do not use on the head or any hard or bony area of the body. Bruising can occur regardless of control setting or pressure applied. Check treatment areas frequently and immediately stop at the first sign of pain or discomfort. NEVER use directly on swollen or inflamed areas or skin eruptions. Batteries (battery pack or batteries installed) shall not be exposed to excessive heat such as sunshine, fire or the like. Danger of explosion if battery is incorrectly replaced. Replace only with the same or equivalent type.

3 YEAR GUARANTEE

FKA Brands Ltd guarantees this product from defect in material and workmanship for a period of 3 years from the date of purchase, except as noted below. This FKA Brands Ltd product guarantee does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorised accessory; alteration to the product; or any other conditions whatsoever that are beyond the control of FKA Brands Ltd. This guarantee is effective only if the product is purchased and operated in the UK / EU. A product that requires modification or adaptation to enable it to operate in any country other than the country for which it was designed, manufactured, approved and / or authorised, or repair of products damaged by these modifications is not covered under this guarantee. FKA Brands Ltd shall not be responsible for any type of incidental, consequential or special damages. To obtain guarantee service on your product, return the product post-paid to your local service centre along with your dated sales receipt (as proof of purchase). Upon receipt, FKA Brands Ltd will repair or replace, as appropriate, your product and return it to you, post-paid. Guarantee is solely through HoMedics Service Centre. Service of this product by anyone other than the HoMedics Service Centre voids the guarantee. This guarantee does not affect your statutory rights. For your local HoMedics Service Centre, go to www.homedics.co.uk/servicecentres

Battery directive

 This symbol indicates that batteries must not be disposed of in the domestic waste as they contain substances which can be damaging to the environment and health. Please dispose of batteries in designated collection points.

WEEE explanation

 This marking indicates that this product should not be disposed with other household wastes throughout the EU.

To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

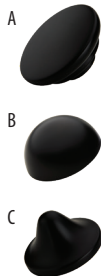
Fig. 1

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FKA Brands Ltd, Somerhill Business Park, Tonbridge, Kent TN11 0GP, UK
FKA Brands Ltd, 29 Earlsfort Terrace, Dublin D02AY28, Ireland