
## SHIATSU DELUXE NECK & SHOULDER MASSAGER WITH HEAT

---

instruction manual

## PRODUCT FEATURES: Convenient and Portable. Perfect for use at home or in the office.
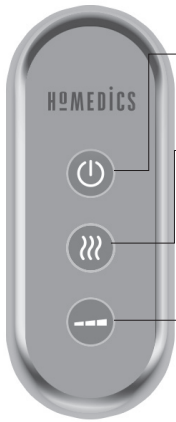
### POWER BUTTON

To turn on the massage function, press the power button once. Press a second time to switch off.

### HEAT BUTTON

For soothing heat, when massage is on, simply press the heat button and the massage nodes on the neck unit will illuminate. To deselect, press the button again and the corresponding massage nodes will stop glowing.

### SPEED CONTROL

The massager will start up at its highest speed setting. Press the speed control button during use to change the speed between the 3 settings - high, med and low.

The actual product materials and fabrics may vary from the picture shown.

Fig. 1

Fig. 2

## INSTRUCTIONS FOR USE:

1. Plug adaptor into a 220-240V AC outlet and plug the receptacle end into the input jack, located along the handle of the neck massager.
2. Place massager around neck (Fig. 1) or area of back you wish to massage. Convenient straps help assist in easy positioning and intensity control (Fig.2).

### **Note: Do not use this product for more than 15 mins at a time.**

There is a 15 minute auto shut-off on this appliance for your safety. This feature should not be considered as a substitute for 'OFF'. Always remember to switch the appliance off when not in use.

## MAINTENANCE:

### **Cleaning**

Unplug the appliance and allow to cool before cleaning. Clean only with a soft, slightly damp sponge.

- Never allow water or any other liquids to come into contact with the appliance.
- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, glass / furniture polish, paint thinner etc to clean.

### **Storage**

Place the appliance in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.

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Distributed in the UK by: FKA Brands Ltd, Somerhill Business Park, Tonbridge, Kent TN11 0GP, UK
EU Importer: FKA Brands Ltd, 29 Earlsfort Terrace, Dublin 2, Ireland