HOMEDICS®

LEADERS IN HEALTH & WELLBEING

MASSEUR 3 EN 1

SHIATSU & THAÏ

SHIATSU AND THAI 3 IN 1 MASSAGER

DOS & EPAULES

BACK & SHOULDER

manuel d'instructions

ASGM-7000ATNDY

GARANTIE 3 ANS

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Unique Gel nodes provide a smooth, natural touch massage.

Demo Button

Activates a brief demonstration of product features.

Up/Down Shoulder adjust $\Delta \nabla$

First ensure neck support pillow is lifted and you are seated before adjusting the height. To adjust height of shoulder unit press and hold Up/Down button (the nodes will partially inflate – hold until the nodes reach your shoulders then release the button – the nodes deflate before massage starts)

Shoulder Air Massage $\curvearrowright \curvearrowleft$

The air pockets will gently begin to inflate giving a relaxing Thai compression massage.

Lumbar Air Massage $\curvearrowright \curvearrowleft$

The air pockets will gently begin to inflate giving a relaxing Thai compression massage.

Intensity $\text{---} \text{---} \text{---}$

3 levels of intensity shown by LED indicators – continue pressing to set the intensity.

Gentle Heat SSS

For gentle soothing heat on Gel Shiatsu back or Thai air shoulder massage, press the gentle heat button. To deselect, press the button again.

Power Button

To start a massage, press the 'POWER' button and the LED will illuminate to confirm your selection. To turn off, simply press the button again. The LED will blink while the massage mechanism returns to its starting position.

Shiatsu Gel Massage

A deeply relaxing and rejuvenating experience invented in Japan using rhythmic finger pressure to stimulate 'chi' (vital energy), relieving stress and tension.

Shiatsu Massager

Full back gel massage

LED indicates selection.

Mid back gel massage

LED indicates selection.

Lower back gel massage

LED indicates selection.

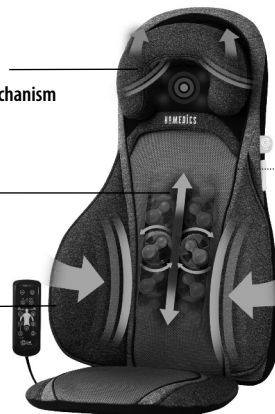
Press again: Stops shiatsu massage.

Thai Air Compression shoulder massage mechanism (lift before use)

Moving Gel Shiatsu massage mechanism

Thai Air Compression lumbar massage mechanism

Aroma tab

The actual product materials and fabrics may vary from the picture shown.

NOTE: There is a 15-minute auto shut off on the unit for your safety.

NOTE: The shiatsu massage mechanism on the back always "parks" or finishes in its lowest position. It will continue to this position after the power has been turned off. If the electrical supply is interrupted, when power is restored, the mechanism will "park" or move to its lowest position. The Thai Air shoulder mechanism will "park" in its uppermost position.

INSTRUCTIONS FOR USE:

IMPORTANT - Preparing the appliance for use

The Shiatsu Massager comes with a screw located on the back of the unit to protect the massage mechanism during shipment, which must be removed before the first use to allow the shiatsu feature to function. Use the provided Allen key. Then, properly dispose of the screw.

Failure to remove screw may cause permanent damage to the appliance.

WARNING Error: all remote control LEDS flash for a few seconds.

Cause: massager overloaded.

Solution: ensure the allen key bolt/screw has been removed from back of the massage cushion (using the allen key included). Should the condition continue after this please contact your HoMedics service centre.

Strapping System

The massager incorporates a unique strapping system which allows you to fasten it to almost any chair. Just slip the straps over the seat or chair and adjust fastening straps to secure. Your massager won't slip or slide away.

NOTE: Do not use this product for more than 15 mins at a time.

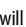
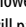
NOTE: complete steps 1 and 2 before switching the appliance on at the mains.

1. Attach the Massage seat to almost any chair using the integrated strap located at the back of the seat. Ensure it is held firmly in place by adjusting the strap as necessary.
2. Connect the power supply lead into a 220-240V AC mains outlet and switch on.
3. Once seated, use the remote control to operate the appliance. Press the  button.

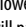
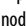
Before use lift the Thai air compression shoulder massage mechanism so you can rest your head comfortably.

Demo mode

- Press demo button
- Up/Down   arrows will flash for shoulders
- Press down  arrow to lower the shoulder mechanism on to your shoulders (the nodes will partially inflate before the mechanism starts to move, hold until the nodes reach your shoulders then release the button - the nodes will then deflate before massage starts)
- Demo mode will then cycle all massage modes

Air Shoulder massage mode

- Use Up/Down arrows  
- Press down arrow  to lower the shoulder mechanism on to your shoulders (the nodes will partially inflate before the mechanism starts to move, hold until the nodes reach your shoulders then release the button - the nodes will then deflate)
- Press the shoulder massage button  to activate Thai air compression massage - intensity can be set using the +/- button 

(Note: intensity for both shoulder and lumbar are controlled simultaneously)

Gentle heat is controlled simultaneously on gel shiatsu back massage and Thai air compression shoulder massage

4. Select lumbar air massage  and/or shiatsu gel massage  as required.
5. Press the  button to stop the massage.

Using essential oils with the aroma tab

Gently pour a few drops of oil onto the aroma tab, releasing the scent in the air. The more drops, the stronger the scent.

MAINTENANCE:

Cleaning

Unplug the appliance and allow to cool before cleaning. Clean only with a soft, slightly damp sponge.

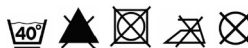
- Never allow water or any other liquids to come into contact with the appliance.
- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, glass / furniture polish, paint thinner etc to clean.

Storage

Place the appliance in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.

Back flap

This massage cushion incorporates a removable back flap in suede finish fabric to add comfort to the massage experience - this is used to soften the massage if too intense. The flap is easily removed as secured by hook and eye fastening. To clean back flap only simply follow the wash instructions:

FKA Brands Ltd

Manufacturer & UK Importer: FKA Brands Ltd, Somerhill Business Park, Tonbridge, Kent TN11 0GP, UK

EU Importer: FKA Brands Ltd, 29 Earlsfort Terrace, Dublin 2, Ireland

Customer Support: +44(0) 1732 378557 | support@homedics.co.uk

IB-ASGM7000ATNDY-1022-02