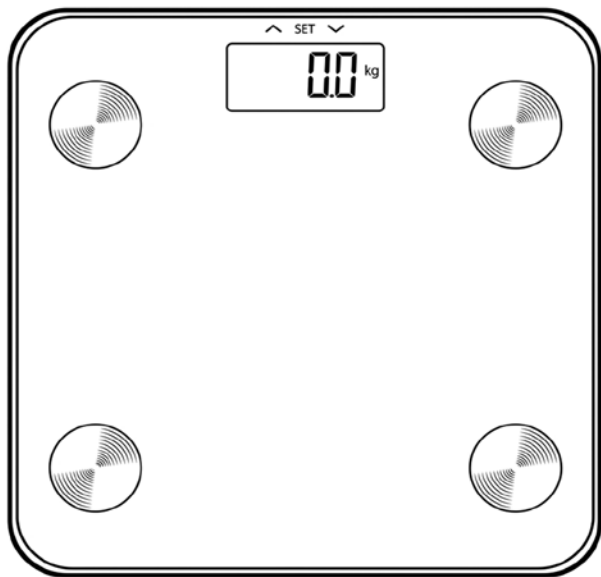
User manual

Scale

To register for your extended guarantee, scan the QR code or visit guarantee.upplc.com/salter and register your product within 30 days of purchase.
For any further queries, visit www.salter.com.

SALTER
1760

EN | Please retain instructions for future reference.

SAFETY INSTRUCTIONS

BATTERY WARNING: KEEP OUT OF THE REACH OF CHILDREN.

- This product is not a toy.
- Children must be supervised to ensure that they do not handle or play with the product.
- This product contains no user serviceable parts. Only a qualified electrician should attempt repairs.
- Do not insert fingers into any openings or inside the product.
- Do not place, store or use the product near direct sunlight or direct heat sources or in a damp/high humidity environment.
- This product is a high precision device; handle it with care.
- This product is intended for domestic use only.

CAUTION: The user should not use the analyser function if they have a pacemaker or other medical device fitted. If in doubt, consult a doctor.

BATTERY SAFETY

- This product is supplied with 2 x CR2032 batteries. Only batteries of the same or equivalent type as those recommended should be used with this product.
- Batteries are to be inserted with the correct polarity.
- Exhausted batteries should be removed from the product to avoid leakage.

	BATTERY WARNING THIS PRODUCT CONTAINS A BUTTON/COIN BATTERY. DO NOT SWALLOW. SWALLOWING CAN CAUSE SEVERE OR FATAL INJURIES WITHIN 2 HOURS. KEEP BATTERIES OUT OF REACH OF CHILDREN (WHETHER NEW OR USED). SEEK IMMEDIATE MEDICAL ATTENTION IF A BATTERY IS SUSPECTED TO BE SWALLOWED OR INSERTED INSIDE ANY PART OF THE BODY.	
--	---	--

CAUTION: Vigilance must be maintained with flat or spare button or coin batteries in the home and the products that contain them.

WARNING: If the battery compartment does not close securely, stop using the product and keep it away from children.

There are no obvious specific symptoms when a button or coin battery is stuck in a child's oesophagus (food pipe). However, the child might exhibit the following symptoms:

- Coughing, gagging or drooling
 - Appear to have a stomach upset or virus
 - Vomiting
 - Pointing to their throat or stomach
 - Pain in their abdomen, chest or throat
 - Tiredness and lethargy
 - Be quieter or more clingy than usual
 - Loss or reduction of appetite
 - Refusal or inability to eat solid foods
- These sorts of symptoms may fluctuate.

A specific symptom of button or coin battery ingestion is vomiting fresh (bright red) blood. Seek immediate medical attention if this occurs.

Care and maintenance

Wipe the exterior of the scale with a soft, damp cloth and allow it to dry thoroughly. Never immerse the scale in water or any other liquid. Do not use harsh or abrasive cleaning detergents or scourers to clean the scale, as this could cause damage.

Before first use

Loosen the battery safety screw and open the battery compartment on the underside of the scale. Remove the plastic isolator tab from the battery compartment (if present) prior to use. Check that the batteries are inserted following the correct polarity. Close the battery compartment and tighten the battery compartment screw. Position the scale on a firm flat surface.

Switching on the scale

Press the centre of the platform to turn on the scale.

Changing the unit

Select kg, lb or st weight mode using the button on the underside of the scale. The scale will then turn off.

Setting the scale to zero

STEP 1: Press the centre of the platform to turn on the scale. Make sure no weight is on the scale. '0' will be displayed.

STEP 2: The scale will switch off after approx. 10 seconds and is now ready for use.

NOTE: The zero-setting process must be repeated if the scale is moved. At all other times step straight on the scale.

Using the scale to take weight readings

STEP 1: Step on to the platform and stand still.

STEP 2: After approx. 2–3 seconds the weight reading will be displayed.

STEP 3: Step off the platform; the weight reading will be displayed for several seconds before the scale turns off.

NOTE: It is recommended that readings are taken at the same time of day, preferably early evening and before a meal for the most consistent results.

Setting user memory

STEP 1: Press 'SET' to turn the scale on. The button is located at the top of the scale behind the glass platform.

STEP 2: Press 'SET' again to enter data entry mode.

STEP 3: Whilst the user number is flashing, select the desired user number by pressing or . Press 'SET' to confirm the selection.

STEP 4: or will flash. Press or to select or . Press 'SET' to confirm the selection.

STEP 5: The height display will flash. Press or to select the correct height, then press 'SET' to confirm the selection.

STEP 6: The age display will flash. Press or to select the correct age, then press 'SET' to confirm the selection.

STEP 7: The display will show the user number and '0.0/ZERO' for approx. 10 seconds, then switch off. The memory is now set.

STEP 8: Repeat the above process to add a second user (max. 9 users), or to change user details.

Using the scale

This scale uses BIA (Bio Impedance Analysis) technology which passes a small electrical impulse through the body to distinguish fat from lean tissue. The electrical impulse cannot be felt and is perfectly safe. Contact with the body is made via stainless steel pads beneath the glass platform.

STEP 1: Position the scale on a firm, flat surface. Press 'SET' to turn the scale on.

STEP 2: Select user number by using or .

STEP 3: The LCD display will show the saved data for that user, followed by a zero reading.

STEP 4: With bare feet, step onto the platform and stand still. Make sure feet are evenly spaced on the platform.

STEP 5: After approx. 2–3 seconds, the weight reading will be displayed.

STEP 6: Remain standing on the scale while the other readings are taken.

STEP 7: The weight will be displayed followed by body fat percentage, body water percentage, bone mass, muscle mass percentage, kcal (kilocalories) and BMI (body mass index).

NOTE: During pregnancy a woman's body composition changes considerably in order to support the developing child. Under these circumstances, body fat percentage readings could be inaccurate and misleading.

Pregnant women should therefore only use the weight function. For more information, visit salter.com/understanding-your-measurements.

Warning indicators

Lo – Replace battery.

Err – The scale weight exceeds maximum capacity.

Err2 – Out of range, incorrect operation or poor feet contact.

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Salter Guarantee

UK CUSTOMERS ONLY

Subject to the care and use being followed, the Salter guarantee promises to replace, substitute or refund any product that fails due to faulty workmanship or materials. This guarantee does not cover any damage caused by accident or misuse and is in addition to your statutory rights. Should any defects due to faulty workmanship or materials appear within 12 months from the date of purchase, please return the goods with proof of purchase to the retailer.

To be eligible for the extended guarantee, go to guarantee.upplc.com/salter and register your product within 30 days of purchase. For any further queries, contact us at support@salterhousewares.co.uk.

EU CUSTOMERS ONLY

EN | To be eligible for the extended guarantee, go to guarantee.upplc.com/salter and register your product within 30 days of purchase. For any further queries, visit www.salter.com.

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EN | Disposal of Waste Batteries and Electrical and Electronic Equipment

This product is labelled with this crossed out wheel bin symbol in accordance with European Directive 2012/19/EC to indicate that it must not be disposed of with your other household waste. Due to the presence of hazardous substances, mixtures or components, electrical and electronic devices that are not

subject to selective sorting are potentially dangerous to the environment and human health. Please check your local city office or waste disposal service for the return and recycling of this product.

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©Salter trademark. All rights reserved.
Ultimate Products UK Ltd.,
Ultimate Products Europe Ltd.,