
User manual

Scale

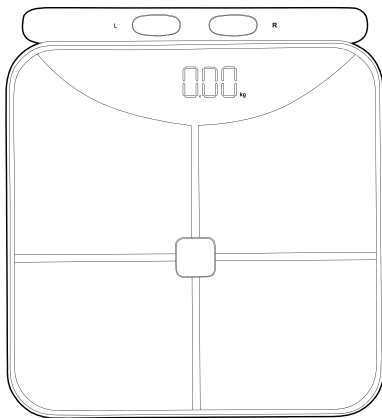For more information, scan the QR code or visit our website: www.salter.com.

Fig. 1

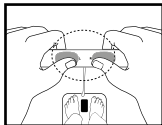
Fig. 2

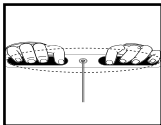
Fig. 3

Fig. 4

To register for your extended guarantee, scan the QR code or visit guarantee.upplc.com/salter and register your product within 30 days of purchase. For any further queries, visit www.salter.com.

SALTER
1760

SAFETY INSTRUCTIONS

BATTERY WARNING: KEEP OUT OF THE REACH OF CHILDREN.

- This product is not a toy.
- Children must be supervised to ensure they do not handle or play with the product.
- This product contains no user serviceable parts. Only a qualified electrician should attempt repairs.
- Do not insert fingers into any openings or inside the product.
- Do not place, store or use the product near direct sunlight or direct heat sources or in a damp/high humidity environment.
- This product is a high precision device; handle it with care.
- This product is intended for domestic use only.

CAUTION: The user should not use the analyser function if they have a pacemaker or other medical device fitted. If in doubt, consult a doctor.

BATTERY SAFETY

- This product is supplied with 4 x AA batteries. Only batteries of the same or equivalent type as those recommended should be used with this product.
- Batteries are to be inserted with the correct polarity.
- Exhausted batteries should be removed from the product to avoid leakage.

Care and maintenance

Wipe the exterior of the scale with a soft, damp cloth and allow it to dry thoroughly.

Never immerse the scale in water or any other liquid. Do not use harsh or abrasive cleaning detergents or scourers to clean the scale, as this could cause damage.

Before first use

Remove the plastic isolator tab from the battery compartment (if present) prior to use.

Make sure the batteries are inserted following the correct polarity.

Position the scale on a firm flat surface.

Switching on the scale

Press the centre of the platform and make sure no weight is on the scale. '0.00' will be displayed.

Converting weight units

Select kg, st or lbs weight mode in the app. The scale will then turn off and is now ready for use.

NOTE: The setup process must be repeated if the scale is moved. At all other times, step straight onto the scale.

Pairing a mobile device

STEP 1: Download and install the Salter Balance app on a mobile device.

STEP 2: Ensure pairing capabilities are enabled on the mobile device.

STEP 3: Open the Salter Balance app and follow the on-screen instructions to register an account or log in to an existing account.

STEP 4: Once an account is set up, navigate to 'Device' in the app.

STEP 5: Select 'Click here to add a device', then select model 'SA00755'. Once paired, a pop up saying 'Binding successful' will be displayed.

Using the scale

This scale uses BIA (Bio Impedance Analysis) technology which passes a small electrical impulse through the body to distinguish fat from lean tissue. The electrical impulse cannot be felt and is perfectly safe. Contact with the body is made via ITO coating on the glass platform and the stainless steel pads.

STEP 1: Open the Salter Balance app once paired with the scale.

STEP 2: Step on the scale with bare feet. Ensure the soles of the feet are in full contact with the top and bottom squares on the scale; these are the electrodes. Ensure there is no contact between the left and right lower limbs.

STEP 3: Lift the handle from the scale, ensuring fingers are placed over the metal electrodes on the back and the thumbs are positioned in the thumb rests on the front. See 'Fig. 1' and 'Fig. 2'.

STEP 4: Keep arms straight at a 30–45° angle. Arms should be straight and not too close to the body. See 'Fig. 3' and 'Fig. 4'.

STEP 5: Maintain this posture until the app prompts that the reading is complete. Once the data is displayed, step off the scale.

STEP 6: When the measurement is complete, the measured weight, body fat, BMI, body water, muscle mass, bone mass and visceral fat will be displayed. Swipe down to view the values of all the listed metrics. Tap 'Report' to view and share the full body analysis report.

STEP 7: To view trends, click 'Trends'.

NOTE: During pregnancy a woman's body composition changes considerably in order to support the developing child. Under these circumstances, body fat percentage readings could be inaccurate and misleading. Pregnant women should therefore only use the weight function. For more information, visit [salter](https://www.salter.co.uk).

com/understanding-your-measurements.
To receive data, the Salter Balance app must be open. The Salter Balance app will not receive data if it is closed or not paired with the scale.

Warning indicators

Lo – Replace the batteries.

Err – Weight exceeds maximum capacity. Unstable weight. Stand still.

C – An error has occurred during measurement. Step off the scale and restart the process.

Specifications

Product code: SA00755

Display: LED

Battery: 4 x AA batteries (included)

Operating frequency: 2.4 GHz

Max RF output: <20 dBm

Units: kg/st/lb

Declaration of conformity

Ultimate Products UK Ltd. hereby declares that this product is in compliance with the essential requirements and other provisions of directive 2014/53/EU.

The declaration of conformity may be consulted at

Ultimate Products UK Ltd., UK. Victoria Street, Manchester OL9 0DD/ Ultimate Products Europe Ltd., 19 Baggott Street Lower, Dublin, D02 X658 or online at **salter.com/declaration-letters/**. Model No. SA00755

For more information, scan the QR code or visit our website: **www.salter.com**.

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Salter Guarantee

UK CUSTOMERS ONLY

Subject to the care and use being followed, the Salter guarantee promises to replace, substitute or refund any product that fails due to faulty workmanship or materials. This guarantee does not cover any damage caused by accident or misuse and is in addition to your statutory rights. Should any defects due to faulty workmanship or materials appear within 12 months from the date of purchase, please return the goods with proof of purchase to the retailer. To be eligible for the extended guarantee, go to **guarantee.upplc.com/salter** and register your product within 30 days of purchase. For any further queries, contact us at **support@salterhousewares.co.uk**.

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**EN | Disposal of Waste Batteries and Electrical and Electronic Equipment**

This product is labelled with this crossed out wheel bin symbol in accordance with European Directive 2012/19/EC to indicate that it must not be disposed of with your other household waste. Due to the presence of hazardous substances, mixtures or components, electrical and electronic devices that are not subject to selective sorting are potentially dangerous to the environment and human health. Please check your local city office or waste disposal service for the return and recycling of this product.

health. Please check your local city office or waste disposal service for the return and recycling of this product.

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©Salter trademark. All rights reserved.