

# Simple Scan Forehead Thermometer/ Touch Free Forehead Thermometer

User's Guide TH106

Please be sure to read all warnings, instructions and notations before using this product. Please keep these instructions as they contain important information.

## Basics of Your Thermometer
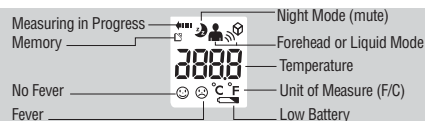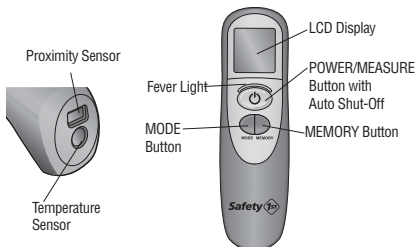
## Care and Maintenance

### Cleaning Temperature Sensor Lens:

- Clean very gently with a soft, alcohol-moistened cloth.
- Do not use an excessive amount of alcohol.
- Do not pour alcohol directly on sensor.

### Cleaning Thermometer Body:

Use a soft, dry cloth to clean thermometer body.

- **DO** remove battery if stored for a long period of time.
- **DO** store thermometer in a dry place at room temperature.
- **DO** take special care not to scratch the sensor or display window.
- **DO NOT** use abrasive cleaning agents, thinners or benzene for cleaning.
- **DO NOT** expose thermometer to extreme temperatures, humidity or direct sunlight.
- **DO NOT** drop thermometer.
- **DO NOT** submerge any part of the product in water.

**PLEASE NOTE:** The sensor lens is the most delicate part of the thermometer. To ensure accurate readings it must be clean and intact.

## Normal Body Temperature

Temperature readings are different from person to person and also vary depending upon the body location from which the temperature is measured. Forehead temperatures tend to be comparable to oral temperatures and approximately 1" lower than rectal temperatures for the same individual. Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.

"Normal" temperature refers to each individual's average temperature while healthy. You may want to measure and record your child's temperature a few times when he or she is well, so that you will know what is normal for your own child. On average, temperature ranges are as follows:

TEMPERATURE METHOD	NORMAL TEMPERATURE	FEVER RANGE
Forehead	36.1°C to 37.9°C	38°C +

normal forehead temp range

## 4 Measuring Modes

Press MODE button to switch between 4 different measuring modes (in order): Forehead , Liquid , Forehead/Night , and Liquid/Night Mode.

NOTE: The beep sounds will be off (mute) when your choice is Forehead/Night Mode, and Liquid/Night Mode.

# WARNING

• This is **NOT** an ear thermometer. Take temperature measurements from the center of forehead. • Use caution when taking temperature readings of bath water and liquids. This product is not intended to touch the person or item whose temperature is being measured. • **ALWAYS** supervise your child when using the thermometer. **NEVER** allow your child to use it alone. • **NEVER** allow children to play with this thermometer. It is not a toy. This is a sensitive electronic device and may be damaged if not used properly under parental supervision. • **DO NOT** expose the thermometer to extreme temperatures, humidity, direct sunlight or shock, and/or dust. • Keep battery out of child's reach. Swallowing a battery could be fatal. Battery should not be charged or placed into extreme heat as they may explode. • Used batteries should be disposed of properly. **Other Battery Information:** We recommend alkaline batteries for longer battery life. **DO NOT** mix old and new batteries. **DO NOT** mix alkaline, standard (carbon zinc) or rechargeable (nickel-cadmium) batteries. Non-rechargeable batteries are not to be recharged. Rechargeable batteries are to be removed from the unit before being charged. Rechargeable batteries are only to be charged under adult supervision. Only batteries of the same or equivalent type as recommended are to be used. Batteries are to be inserted with the correct polarity. The supply terminals are not to be short-circuited. Check that all contact surfaces are clean and bright before installing batteries. **DO NOT** disassemble for any reason other than inserting or replacing batteries.

## How to Measure Temperature


Please read this section carefully and all the way through before using the thermometer. Ensure that thermometer is at room temperature, waiting 15 minutes if necessary.

**Turn On:** Press and release the POWER/MEASURE button . All symbols are displayed (Figure A).

**Verify/Select Mode:** The default setting is "Forehead Mode". Verify that the thermometer is in the correct mode — "Forehead Mode" (Figure B) or "Liquid Mode" (Figure C)— for the temperature you will measure. Press MODE button to change current setting.

## Forehead Temperature

1. Position sensor approximately 1" from center of forehead (between eyebrows) (Figure D).
  2. Press the POWER/MEASURE button. The symbol, measurement in progress symbol and "00" should be displayed (Figure E).
  3. Move thermometer closer or farther from the forehead until you hear a short beep and see the screen light up. This means the temperature reading is complete.
  4. If the temperature is in normal range a NO FEVER icon will appear next to the reading. If elevated temperature is detected, a FEVER icon will be displayed and the fever light above the POWER/MEASURE button will illuminate orange. **NOTE: High, prolonged fever requires medical attention, especially for young children. Please contact your physician.**
  5. If the thermometer is too far from forehead, a dash-icon will flash on screen accompanied by a beep sound until the correct distance is achieved. This will last for 20 seconds and if after 20 seconds the right distance is not achieved, the display will turn to stand-by mode (Figure F).
- NOTE: Additional temperature measurements can be taken continuously.
6. Press the POWER/MEASURE button to turn thermometer off. If left on, it will automatically shut off in about 30 seconds.
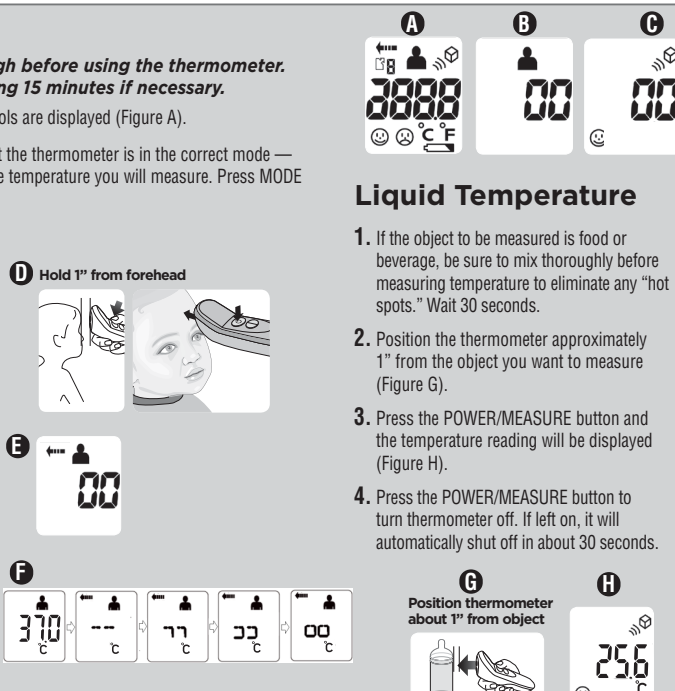
## How to Install/Replace the Batteries


This thermometer comes with the batteries outside of the unit. If this is your first time using this unit, start with step 1 below.

**Low Battery Warning:**

- When the symbol and the icon appear on the display screen, the batteries need to be replaced before using thermometer again.

**To Install Batteries:**

- 1 Slide cover downward to remove cover.
- 2 Remove old batteries.
- 3 Insert two AAA (or equivalent) batteries, correctly positioned with positive "+" and negative "-" as shown in illustration.
- 4 Slide the battery cover onto the thermometer.
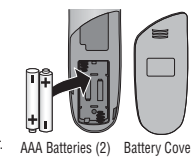
## Memory

The thermometer memory function stores the last ten measurements. Press and release the MEMORY button to see previous readings stored in memory. Each press displays a previous reading with the Forehead or Liquid symbol.

NOTE: Press and hold the MEMORY button for more than 3 seconds to delete all stored readings. Removing the batteries will clear all stored memory.

## Switch Between Fahrenheit & Celsius

While the thermometer is ON, press and hold the MODE and MEMORY buttons for about 1 second until you see the °F or °C. Release buttons.

**Press 2 buttons**

## Technical Specifications

ASTM laboratory accuracy requirements in the display range of 36 to 39 °C (96.8 to 102.2 °F) for IR thermometers is +/- 0.2 °C (+/- 0.4 °F), whereas for electronic thermometers, the requirement per ASTM Standard 1965-98 is +/- 0.1 °C (+/- 0.2 °F).

**Measurement Temperature Range:**  
Object: 32°F to 212°F Human Body: 93.2°F to 109.4°F

**Measurement Accuracy:**  
+/- 0.4°F Temps 96.8°F to 102.2°F  
+/- 0.5°F All Other Temps

**Operation Environment:**  
59°F to 95°F  
Humidity up to 95% (non condensing)

**Storage Environment:**  
-4°F to 122°F  
Humidity up to 85% (non condensing)

**Test Protocol:**

- 1 Place the IR thermometer into water bath chamber for at least 2 hours.
- 2 After 2 hours, press MEASURE button without releasing key until the decimal point disappears.
- 3 Put the thermometer into the fixture of black body.
- 4 Press and release MEASURE button, the device will beep once and measurement will be shown on the display.

## One (1) Year Limited Warranty

Dorel Juvenile Group, Inc. warrants to the original purchaser that this product (Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc., will repair or replace the product, at our option, free of charge. Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on the front of this document and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse. This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or consequential damages.

**SOME PROVINCES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.**

**THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS THAT MAY VARY FROM PROVINCE TO PROVINCE.**

Manufactured by/Fabriqué par: Dorel Juvenile Group, Inc. A division of/Une division de Dorel Industries, 2855 Argentinia Road, Unit 4 Mississauga, ON L5N 8G6 Made in CHINA. Fabriqué en CHINE. (800) 544-1108 www.safety1st.ca Styles and colours may vary. Les styles et couleurs peuvent varier. Distributed by/ Distribué par: Dorel Juvenile Group, Inc., 2525 State St., Columbus, IN 47201-7494.

Dorel Juvenile Canada, 2855 Argentinia Road, Unit 4, Mississauga, ON, L5N 8G6 01/06/21 4358-8097B

Display Meaning		Recommended Action	Troubleshooting		Display Recommended Action
	Measured temperature is below this thermometer's range. See Technical Specifications section for Measurement Temperature Range.	<ul style="list-style-type: none"> <li>• Ensure thermometer is being positioned correctly (See How to Measure Temperature).</li> <li>• Ensure  icon does not appear. If it does, replace batteries.</li> <li>• Consistent  readings may indicate that the sensor is dirty (see Care and Maintenance).</li> </ul>		Batteries are low. Replace them. (See Battery Information section).	
	Measured temperature is above this thermometer's range. See Technical Specifications section for Measurement Temperature Range.	<ul style="list-style-type: none"> <li>• Ensure thermometer is being positioned correctly (See How to Measure Temperature).</li> <li>• Allow thermometer to reach room temperature, then re-take temperature measurement.</li> </ul>		<ul style="list-style-type: none"> <li>• For best accuracy, wait one minute between temperature measurements so that thermometer has time to re-acclimate to the room's temperature.</li> </ul>	
	Thermometer may have been stored in a cold or hot environment. This thermometer cannot function if the environment's temperature is outside the permissible range. See Operation Environment for range.	<ul style="list-style-type: none"> <li>• Allow thermometer to reach room temperature, then try again.</li> </ul>		<ul style="list-style-type: none"> <li>• Do not return this product to the place of purchase. To address any issues please visit us at www.safety1st.ca to chat with us there. You can also text us at (812) 652-2525 or call Consumer Care at (800) 544-1108. Have ready the model number (TH106) and date code (manufacture date located on package).</li> <li>• Install new batteries and re-try function before calling Consumer Care on electronics.</li> </ul>	

Model name: NT11  
[REF] TH106  
Manufacturer: AVITA (Wujiang) Co. LTD.  
Tel: +86-512-8207-8288  
执行标准: GB/T 21417.1-2008  
医疗器械注册证号: 苏械注准20182200808  
医疗器械许可证号: 苏食药监械生产许20040018号  
Add.: No. 858, Jiao Tong Rd., Wujiang Economic Development Zone, Suzhou City, Jiangsu Province, P.R.C.

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