Scan for support  
Scannez pour l'aide  
Escanee para ayuda

**1-812-203-4747**

(US and Canada only, carrier fees may apply)  
(Uniquement pour les États-Unis et le Canada, des frais peuvent s'appliquer)  
(Sólo para EE.UU. y Canadá, pueden aplicarse tarifas de operador)

**1-800-544-1108**


**safety1st.com/connected-support**

**connected@safety1st.com**

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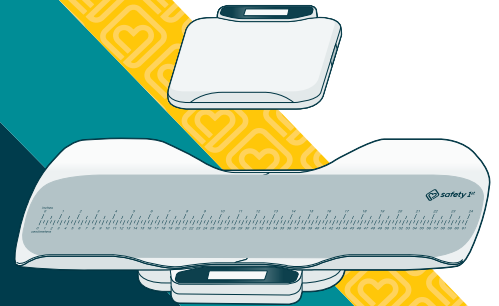## Smart Family Scale Balance intelligente Balanza inteligente

Install the App  
Installez l'application  
Instale la aplicación

Log in or create an account,  
then tap **+** to set up your scale.

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Keep these instructions for future reference.

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**Weight range:**

- For baby tray: 0 to 44 lbs (0–20 kg) or 24 months, whichever comes first
- For scale: 0 to 330 lbs (0–150 kg)

**Measurement precision:**

- +/- 0.71 oz up to 22 lbs (+/- 20 g up to 10 kg)
- +/- 7.06 oz up to 330 lbs (+/- 200 g up to 150 kg)

# ⚠ WARNING

- Always use the scale on a solid, flat surface. DO NOT use on soft ground, for instance: on a blanket, mattress, or carpet.
- Do not knock, impact, or drop the scale.
- This product is a precision instrument, please be careful.
- The platform will be slippery when it is wet. Keep it dry at all times.
- Keep the scale in cool place with dry air.
- If the scale fails to turn on, ensure that batteries are properly installed or whether the battery power is low. Insert or replace batteries as needed.
- Designed for family and non-commercial use.
- Keep still during the course of weighing.
- Adult assembly is required. Keep small parts away from children when assembling. Exercise care when unpacking the product.
- This product is not a toy. DO NOT allow children to play with it.
- Make sure all tray fasteners are fully secured to the scale before using for your baby.
- Clean the scale with a damp cloth and avoid water getting inside the scale. DO NOT use chemical/abrasive cleaning agents.
- DO NOT wash cork mat in washing machine.
- Dispose of batteries safely.
- This product contains magnets. Swallowed magnets can stick together across intestines causing serious infections and death. Seek immediate medical attention if magnets are swallowed.

**Important Battery Information:**

- Keep these instructions as they contain important information.
- We recommend alkaline batteries for longer battery life.

# ⚠ BATTERY WARNING

**CAUTION: TO AVOID BATTERY LEAKAGE:**

- Always follow the instructions carefully. Use only specified batteries. Make sure to insert batteries with the correct polarity.
- Remove exhausted batteries from the product.
- Do not short-circuit the supply terminals.
- Remove batteries when stored for long periods of non-use.
- **RECHARGEABLE BATTERIES:** Rechargeable batteries are to be removed from the product before being charged. Rechargeable batteries are only to be charged under adult supervision. Do not recharge non-rechargeable batteries.

Recycle or dispose of batteries safely according to federal, state, and local laws.


**Disposal:**

At the end of the product life cycle, you should not dispose of this product with normal household waste. Take this product to a collection point for the recycling of electronic equipment in accordance with your local laws. Some of the product materials can be reused if you take them to a recycling point. By reusing some parts or raw materials from used products, you make an important contribution to the protection of the environment. Please contact your recycling center for more information on the collection points in your area.

## Parts and Features

- a Adult scale
- b Baby tray (2 halves)
- c Baby tray support brackets (2)
- d Cork mat
- e Battery compartment
- f Reset/pairing button (under the scale)

- Tools needed: Philips head screwdriver (not included).
- Requires 4 AAA alkaline batteries (included).

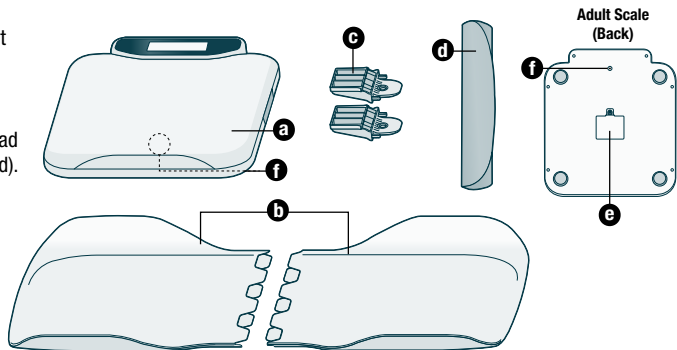
## Display:

Power On

Tare

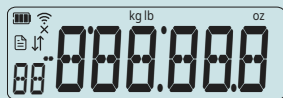
lb kg


Change Units  
(lbs / lbs oz / kg)

## Icons:


Battery Indicator

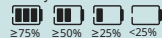Wi-Fi Indicator

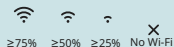Uploading Data Indicator

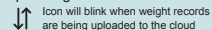Memory Indicator

## Batteries

### Low battery:


Replace batteries when battery LED  indicates low or device doesn't turn on when you press the power button.

### To install or replace batteries:

Use a Phillips head screwdriver to remove the battery door. Remove all batteries and replace with new batteries. Place 4 AAA alkaline batteries inside compartment as shown. Match the (+) and (-) symbols according to the diagram in the battery compartment. Replace the door and secure with the screw.

### Note:

- The unit will not operate if the batteries are installed incorrectly.
- Install new batteries and retry function before contacting Consumer Care.



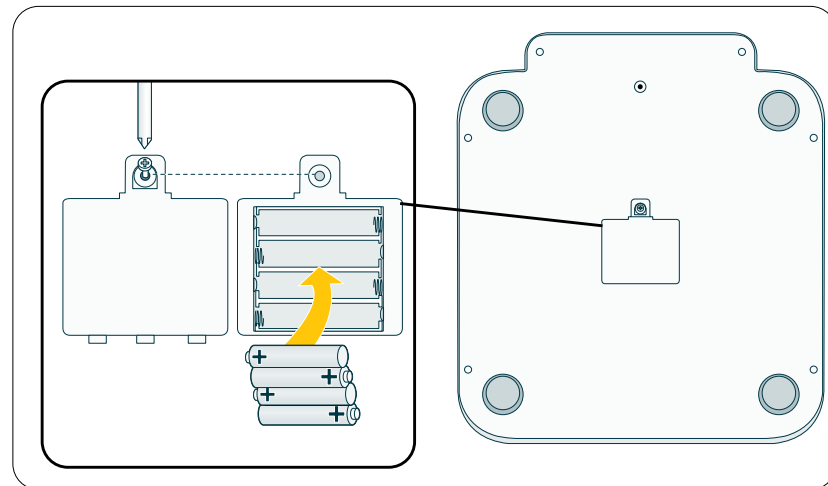
### Assemble Baby Tray

Turn tray halves over and rest on flat/even surface. Orient the brackets as shown. Slide mount brackets into slots on tray until they click in place.

**Tip:** To remove brackets from tray, take the tray off the scale and turn it over. Insert a screwdriver and apply gentle pressure to push the bracket latch into the notch on one side of the tray. While holding that side of the latch (to prevent it from clicking back in), repeat with the other notch until the bracket slides out. Repeat with the other half of the tray.

Locate the opening at each side of the scale. Align the tongue of each mounting bracket with the openings at the sides of the scale. Push tray together until it clicks.

**Note:** The display will rotate 180° when the tray is installed.
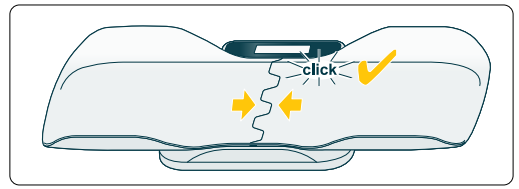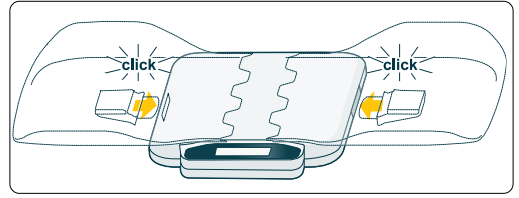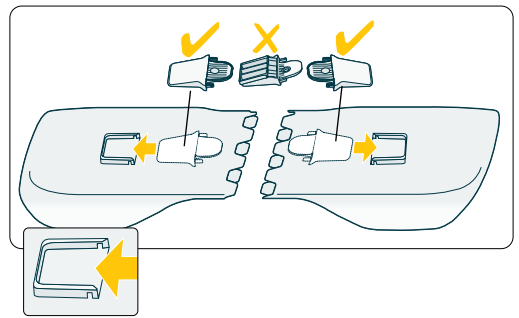
### Setup

1. Download the **Connected Family app**. Ensure that your smartphone is connected to your home's Wi-Fi network.
2. Install batteries (see Page 5).
3. In the app, tap **+**, and select Safety 1st, then tap "Smart Family Scale" to begin pairing your device.
4. Follow step-by-step instructions to complete pairing.
5. Once successfully paired, the Wi-Fi icon will turn on and show signal strength.
6. Tap **?** in the **Connected Family app** to learn more, including how-to videos.

### Wi-Fi Connectivity

For app connectivity, this device requires a 2.4 GHz wireless connection with network speeds of 2 Mbps or greater at the location the device will be installed. Your internet service provider can verify the quality of your internet connection.

If your network equipment supports both 2.4 GHz and 5.0 GHz bands, we recommend that the bands are separated so each have their own Network Name (SSID) and Password (Network Key). Refer to your internet service provider or equipment provider for assistance with this.

### Change Unit of Measurement

**Note:** Factory default is set to pounds. The three available selections include Pounds (decimal), Pounds (with Ounces), and Kilograms.

1. Press the **⏻** On button.
2. Press **⏻** repeatedly until unit of measure desired appears on the display.
3. Wait a moment until your selection blinks. Your selected unit of measurement is ready.

### Reset

Press and hold the reset button for 5 seconds (on the back) until you see the display reset. Unit of measure will be reset to pounds.

### Use Baby Weighing Mode

**Tip:** When weighing baby, position the display nearest to you.

1. Ensure tray is installed and engaged. Scale must be on a flat, stable surface.
2. Place cork mat on tray before turning on.
3. Press the  On button.
4. Place baby on the tray and wait until baby stops moving.  
**Tip:** Weighing baby while they are sleeping may be easier than when they are awake.
5. Once a weight is recorded, the scale will beep and blink briefly before displaying the result for 15 seconds.

### Use Regular Weighing Mode


1. Press the  On button or step on scale.
2. Step on the scale and stand still.
3. Once a weight is recorded, the scale will beep and blink briefly before displaying the result for 15 seconds.

**Note:** Do not stand on the scale with shoes on; wear socks or bare feet.

### Disassemble Baby Tray

When pulling tray halves, pull straight out and then slightly downward. Set tray aside for later use.
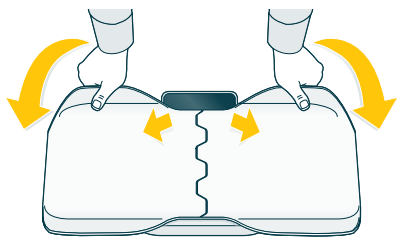
### To Clean

Do not submerge the scale in water. Wipe all parts with damp cloth and dry. Do not wash cork mat in washing machine.

### FCC & IC Statement

This device complies with part 15 of the FCC Rules and Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

### FCC & IC Radiation Exposure Statement:

This equipment complies with FCC and Canada radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located with any other antenna or transmitter.

### RADIO AND TELEVISION INTERFERENCE:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

1. Reorient or relocate the receiving antenna.
2. Increase the separation between the equipment and receiver.
3. Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
4. Consult the dealer or an experienced radio/TV technician for help.

**Note:** Changes or modifications not expressly approved by the manufacturer responsible for compliance could void the user's authority to operate the equipment.

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©2024 Dorel Juvenile Group, Inc. Safety 1st™ and the Safety 1st logo and design are trademarks of Dorel Juvenile Group, Inc. (800) 544-1108 www.safety1st.com. Styles, colors, and labeling may vary. Les styles, les couleurs et l'étiquetage peuvent varier. Los estilos, colores y etiquetado pueden variar. Made in China. Fabriqué en Chine. Hecho en China. Distributed by (distribué par/distribuido por) Dorel Juvenile Group, Inc., 2525 State St., Columbus, IN 47201-7494. Dorel Juvenile Canada, 2233 Argentia Road, 110, Mississauga, ON L5N 2X7. Representante exclusivo en Chile: DOREL JUVENILE CHILE S.A. Av. del Parque 4161, oficina 602 A, Huechuraba, Santiago, Chile - Fono: 56 2 24971500 - www.infanti.cl - www.dorel.cl. Importado y distribuido en Perú por: COMEXA COMERCIALIZADORA EXTRANJERA S.A. SUCURSAL DEL PERÚ - Los Libertadores 455, San Isidro, Lima 27, Lima, Perú RUC 20431062870 - Fono: 4413029 - www.infanti.com.pe. En Centro América, Caribe y Ecuador por: Best Brands Group S.A. - Av. Balboa, Bay Mall Plaza, Planta Baja, local 9, Panamá. - Fono: (507) 300 2884. En Colombia por: Baby Universe S.A.S. - Vía 40 No. 77-29, Barranquilla, Colombia - Fono: (575) 353 1110 Importador: DJGM, S.A. DE C.V. Gabriel Mancera No. 1041 Col. Del Valle, Benito Juárez C.P. 03100 México, Ciudad de México R.F.C. DJG140305GN9 Tel. (55) 67199202 05/24 4358-8489