

medela

disposable inserts

designed for use with the medela waterless milk warmer™

- fast, easy and convenient
- accommodates most human milk storage containers
- accommodates syringes 1 mL - 60 mL



see back panel for instructions

25 - inserts

medela

Waterless Milk Warmer™ and Disposable Inserts

- Warms milk to temperatures consistent with expressed human milk.*
- Eliminates contamination risk of warming milk with tap water.
- Safely thaws milk to a refrigerated temperature for storage up to 24 hours.
- Easy-to-use, easy-to-clean and convenient.
- Inserts contain spills, minimizing clean up time.
- Inserts are single-patient items.



For more information on the Medela Waterless Milk Warmer please visit www.medela.com.

*Temperatures may vary depending on actual container used. Device is optimized for syringes and the Medela Sterile 80 mL Breastmilk Storage Containers. Designed to warm human milk to 90 °C - 38 °C with most commonly used containers found in NICU environments.

medela


Instructions For Use

Review instruction manual for important safeguards prior to first use.

Set-up



Expand liner prior to installing in warmer.

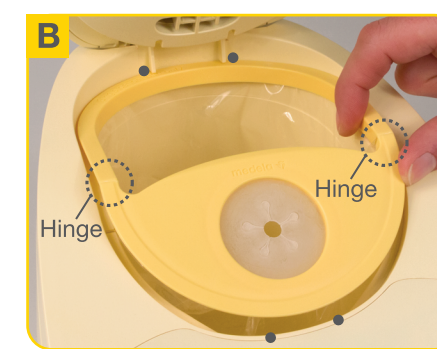

Bend insert at hinge and place into position.



Evenly spread liner within warmer cavity.



Push tray, flat-side down, into bottom of warmer.



Load feed.
NOTE: Always use insert tray to warm container or syringe.

F Press power button to turn on the device.

Continued instructions on side panel

Push Here Then Pull

medela

Program

Select Mode

A Select
WARM
or
THAW
using:
< >
Press enter.

Select Starting Temperature

B Select starting temperature of milk
Refrigerator or Room or Frozen
using:
< >
Press enter.
Note: Thaw Mode automatically defaults to: Frozen

Select Volume


A Select volume of milk to be warmed
000 mL
using:
< >
Press enter.
Note: Hold arrow key to advance more quickly.

Confirm & Begin

B Confirm settings
☒ Mode
☒ Start Temp
☒ Volume
and press
start
button to begin.

Ready

Flashing green light indicates milk is ready.

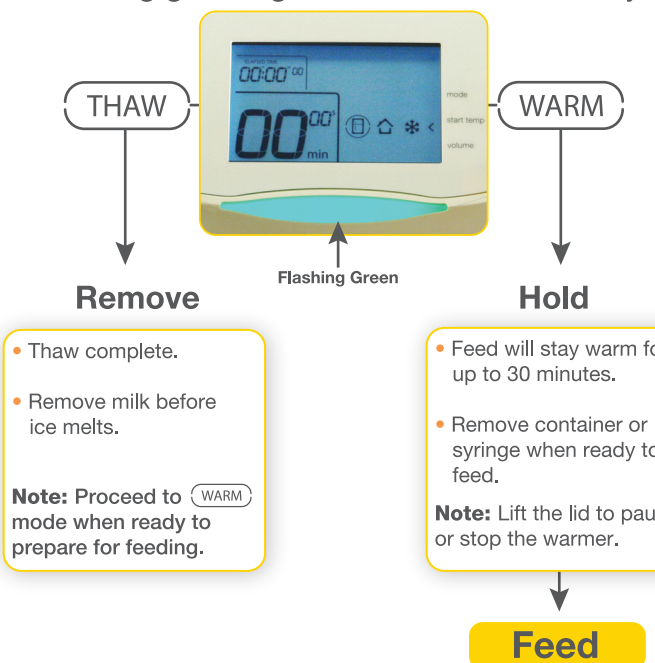

25 - inserts

Content appearance may vary from picture.
Medela, Inc., 1101 Corporate Drive, McHenry, IL 60050
Phone: (800) 435-0316 or (815) 365-1166 Fax: (815) 365-1246
Email: customer.service@medela.com www.medela.com
Packaged and printed in the USA. Manufactured by Medela, Inc.
Medela is a registered trademark of Medela Holding AG.
Waterless Milk Warmer is a trademark of Medela, Inc.
© 2014 Medela, Inc. 1948831 B 0314 All rights reserved.

#87116