



# Oxivir® Plus

30

## Disinfectant Cleaner Concentrate

Dilutable one-step disinfectant cleaner, hydrotherapy tub cleaner and sanitizer, formulated with proprietary hydrogen peroxide technology. Ecologo certified.

### Features & Benefits

- Quick disinfection in just 5 minutes against Human Immunodeficiency Virus (HIV-1) (AIDS Virus), Polio virus, Influenza A Virus, *Pseudomonas aeruginosa*, *Staphylococcus aureus*, *Salmonella enterica*, MRSA and VRE
- 5 minute Bactericidal, Tuberculocidal, Virucidal, Fungicidal with a 30 sec 1:128 non food contact sanitizer claim
- AHP® active breaks down into water and oxygen after use
- No VOCs, no added fragrance and no NPEs
- Delivers fast and effective cleaning performance while being non-irritating to eyes and skin in use dilution

### Applications

- Clean and disinfect non-critical medical devices, equipment and hard, non-porous surfaces





### Use instructions

Oxivir® Plus Disinfectant Cleaner Concentrate is a one-step disinfectant cleaner using patented AHP® technology for use on walls, stainless steel, glazed porcelain, plastic surfaces, hydrotherapy tubs and other hard, nonporous surfaces.

- To disinfect, apply a 1:40 dilution to thoroughly wet surface.
- Let solution remain on surface for a minimum of 5 minutes.
- Air dry, wipe surfaces to dry and remove any residue, or rinse with potable water as necessary.
- Prepare a fresh solution if the solution becomes visibly dirty or diluted.


### USE DIRECTIONS FOR HYDROTHERAPY TUBS AND FOOT BATHS

- To reduce the generation of foam, add an anti-foam agent at the rate of 5% of the amount of this product. Heavily soiled tubs require pre-cleaning prior to disinfection.
- **Initial Purging and Disinfecting** Fill tub with cool water to above the water jet line. Add this product at 1:40.
- Turn on the pump and allow the water to circulate for 5 minutes. Drain then rinse. Broad-Spectrum Sanitizing - Maintenance Fill tub with cool water to above the water jet line. Add this product at 1:128.

- Turn on the pump and allow the water to circulate for 30 seconds. Ensure surfaces have been scrubbed.
- Drain then rinse.

### Dilutable Disinfectant Cleaner

- 1:40 dilution – Bactericidal and virucidal, including antibiotic-resistant bacteria, in 5 minutes.
- 1:128 dilution – Non-food contact sanitizer in 30 seconds or general purpose cleaning.

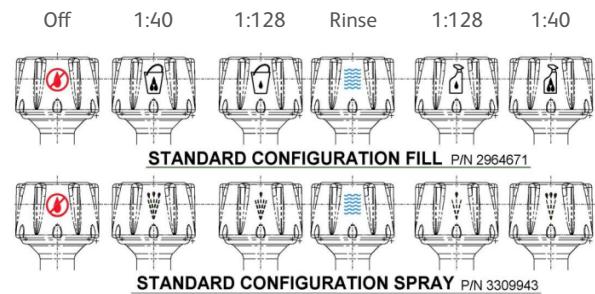

Technical data	Oxivir® Plus
DIN	02403684
Certifications	Ecologo, LONO
Color/Form	Colorless, liquid
pH	0.8 (Concentrate) 1.9 (Use Dilution 1:40)
Scent	Characteristic of ingredients
Shelf Life	2.5 years (Concentrate), 30 days (Use Dilution)

Product	Pack size	Dilution	Product code
Oxivir® Plus	4 x 3.78 L / 1 U.S. Gal. Containers	1:40, 1:128	5919024
Oxivir® Plus	2 x 1.5 L RTD® Units / 50.7 oz.	1:40, 1:128	101109764
Oxivir® Plus	2 x 2.5 L J-Fill® / 84.5 oz.	1:40	101109765

### Safe handling

Please make sure your employees read and understand the product label and Safety Data Sheet before using this product. The label contains directions for use; and both the label and SDS contain hazard warnings, precautionary statements and first aid procedures. SDS are available online at [www.diversey.com](http://www.diversey.com) or by calling 888.352.2249. Improper use or dilution may result in damage to surfaces and may result in health and physical hazards that match those of the concentrate.

Accelerated Hydrogen Peroxide® and Design and Oxivir® are trademarks of Diversey, Inc.