

INVACARE®
PLATINUM® MOBILE
OXYGEN CONCENTRATOR
WITH CONNECTIVITY†

Now with advanced intelligence and five flow settings.



DO MORE WITH PORTABLE OXYGEN



Breathe Smart.
Breathe Easy.


Invacare's Game Changing Oxygen Solution



Piccolo O₂ App from iTunes® (iOS) or Google Play (Android).
Phone and laptop not included.

- ✓ A durable POC with our easiest battery management system and other unique features
- ✓ A user-friendly app[†] that makes time management a breeze
- ✓ A Provider Portal[†] to help you with efficient asset utilization and lower total cost of ownership

All fully integrated to help support outcomes and deliver essential value where it matters the most

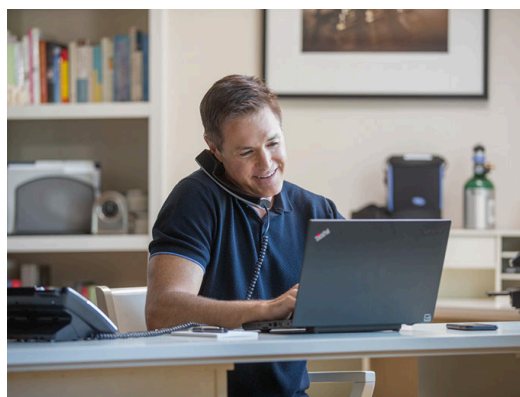

The TOTAL PACKAGE to help you take control of your business performance



Piccolo O₂ App from iTunes® (iOS) or Google Play (Android).
Phone and laptop not included.

What's in it for me?



Asset Utilization and Deployment Management[†]

- Visibility to runtime on deployed units so you can match equipment to the end-user



Lower Total Cost of Ownership[†]

- Portal supports preventative and predictive maintenance
- Helps reduce unplanned truck rolls
- Remote trouble shooting helps increase call center productivity
- End-user self-help app to help support your customers and reduce calls to you



Customer Delight with the Most Unique Features

- Win over your referral sources with a unique, attractive POC designed to support end user compliance
- Delight your customers with a durable POC with our easiest battery management system and other unique features


A World Class Portable Oxygen Concentrator

Lightweight and long battery time are basic* requirements, but your customer deserves more than basic.

- » **3-way convertible bag** (with backpack mode included) for increased mobility
- » Built in protective features for **durability**
- » Easy **battery management** system
 - ✓ Short charge time
 - ✓ Increased battery time, with reduced weight
 - ✓ Easily load battery from top without removing bag
 - ✓ Easily swap batteries without turning off device and interrupting oxygen flow
- » **Smart Piccolo O₂ App[†]**
 - ✓ Shows battery time remaining in hours and minutes not just battery charge level
 - ✓ Premium access to local environmental conditions like humidity, allergies, temperature
 - ✓ Quick access to help, user manuals, lifestyle resources and much more

All this, built into a stylish, easy-to-use unique device.

**See specs page for list of basic features including weight, battery duration, sound pressure level, and more.*

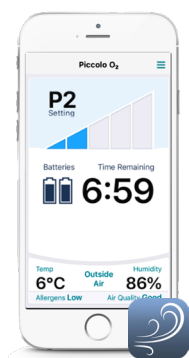

Phone not included.



Desktop battery charger available as a separate purchase (POC1-115).

Designed for Clinical Efficacy

Supplemental oxygen users have varying physiological and respiration profiles. For adequate saturation, they need an oxygen concentrator that delivers the right volume of therapeutic grade oxygen at the right time matching their breathing patterns and with little user effort.

The Invacare Platinum Mobile Oxygen Concentrator with Connectivity† features:

- ✓ **Large bolus volume:** Up to 67ml per bolus*
- ✓ **High oxygen minute volume delivery:** As a function of purity and large bolus size
- ✓ **Quick bolus delivery at inspiration:** Designed to minimize delay and deliver oxygen at inspiration. This covers the most critical part of inhalation and maximizes active gas exchange
- ✓ **Great trigger sensitivity:** Optimal balance between sensing inspiration and avoiding unintended triggering
- ✓ **Invacare Sensi-Pulse™ technology:** Constant oxygen minute volume helps ensure delivery of the highest oxygen purity

**At flow setting 5, 15 BPMs*

All this, built into a stylish easy-to-use unique device, can help make ambulation easier for patients, support compliance, and help reduce costly re-admissions.

DO MORE WITH PORTABLE OXYGEN



PRODUCT SPECIFICATIONS

MODEL #: POC1-100C
HCPCS: E1390/E1392****

Model #	Description
POC1-100C	Concentrator, battery, carrying case, AC adapter and DC adapter
POC1-100CA	Concentrator, dual batteries, carrying case, AC adapter and DC adapter
POC1-110	Supplemental battery - US
POC1-110-C	Supplemental battery - CA
POC1-115	Desktop battery charger- US
POC1-115-C	Desktop battery charger- CA
POC1-130	AC adapter
POC1-140	DC adapter 12V
POC1-150	Carrying case
POC1-CONNECT	Invacare USB Dongle, Platinum Mobile



1. Desktop battery charger
2. 3-way convertible carrying case (included)
3. Easy access HEPA filter
4. Top load & hot swappable battery management system

AVERAGE PULSE VOLUME						
PULSE SETTING	BREATHS PER MINUTE					
	15	20	25	30	35	40
P1	14.7	11.0	8.8	7.3	6.3	5.5
P2	29.3	22.0	17.6	14.7	12.6	11.0
P3	44.0	33.0	26.4	22.0	18.9	16.5
P4	58.7	44.0	35.2	29.3	25.1	22.0
P5	66.7	50.0	40.0	33.3	28.6	25.0

BATTERY DURATION		
PULSE SETTING	ONE BATTERY	TWO BATTERIES
P1	5 HR 5 MIN	10 HR 10 MIN
P2	3 HR 30 MIN	7 HR 0 MIN
P3	2 HR 20 MIN	4 HR 40 MIN
P4-P5	1 HR 45 MIN	3 HR 30 MIN



Dimensions

9.45" H
7.5" W
3.88" D



Weight

5.0 lb.*
15.9 oz. with Battery



Oxygen Flow

1-5 Pulse Settings



Oxygen Purity

87% - 95.6%



Maximum O₂ Capacity

1,000 ml/min
(@ setting 5)



DC Cable Input

12 V
13.5 V Nominal



Limited Warranty

5 Year Concentrator
1 Year Sieve Beds & Battery



AC Adapter Input

100-200 VAC
50/60 Hz
1.5 A - 0.6 A



Battery Duration

One Battery: 3.5 hrs.
Two Batteries: 7 hrs.
(@ pulse setting 2)



Recharge Time**

One Battery: 2 hr 18 min
Two Batteries: 4 hr 36 min



Operating Altitude

Up To 10,000ft. above sea level



Sound Pressure Level

≤40 dB (A) @ Setting 2***



FAA Compliant

Yes



Temperature Tolerance

41° to 104° F
(5° to 40° C)



Water Resistance Rating

IP22 for up to 10 minutes

* Weight with single battery without carrying case.

** Recharge time increases if unit is running while charging.

*** With carrying case.

****All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

Please refer to Owner's Manual for product specifications and product use.



Visit Us Online:

www.Invacare.com/BreatheSmartBreatheEasy
for more information.



† Access to the features and benefits of the Provider Portal and Piccolo O₂ App require Connectivity, which requires purchase of the POC1-CONNECT. Please contact your Invacare Sales Representative or Customer Service for more information.

©2019 Invacare Corporation. All rights reserved. Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Google Play and Android are trademarks of Google, Inc. Windows is a registered trademark of Microsoft Inc. The Bluetooth word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. Apple, iOS, and App Store are registered trademarks of Apple, Inc. Specifications are subject to change without notification. Form No.18-277 rev0119 190112