



INVACARE® MICROAIR® MA1000
THERAPEUTIC SUPPORT SURFACE



Yes, you can.®

microAIR®

Invacare® microAIR® MA1000

Alternating Pressure Low Air Loss Mattress System

**Medium - High Risk
Stage I-IV**

3-in-1 Alternating Pressure System combines
True Low Air Loss and advanced clinical
features, creating the Ultimate Solution for
Patient Comfort and Pressure Relief



Highly vapor permeable and oversized pliable quilted nylon top cover designed to provide low shear, friction and moisture protection

Added convenience of **Integrated side bolsters** that can be inflated and deflated on-demand



Nylon non-slip base cover designed to help reduce shear and friction during bed articulation



Advanced 3-in-1 alternating function inflates and deflates cells in a 3:1 cycle, providing active pressure relief by supporting 2/3 of the body at any given time

Features & Benefits

10 adjustable comfort settings depending on the patient's therapy requirement



Auto Set Mode automatically sets optimal pressure based on patient's weight, height and position for constant pressure monitoring

60 second pulsation feature designed to help improve capillary blood flow by continually readjusting to a level surface pressure.

Auto Fowler Mode safety feature can detect the pressure of patient's buttocks automatically and adjust pressure accordingly

True Low Air Loss Therapy designed to provide extra comfort for moisture and temperature control

Patient Care Mode designed to provide quick maximum inflation within seconds to aid with transfers and nursing procedures

Customizable cycle times to meet patient's varying needs



Foam base designed to help prevent bottoming out in the event of a power failure

Soft nylon urethane backed cells designed to provide easy cleaning and helps minimize the proliferation of bacterial growth

HeelSense™ Technology designed to help provide further therapy and comfort by decreasing pressure in patient's vulnerable heel area

Request a Quote or Place an Order Today!



**Contact your Inside
Sales Representative**

or



**Contact Customer Service
at 800.333.6900**

Invacare Corporation
www.invacare.com

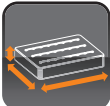





USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Invacare microAIR MA1000 Alternating Pressure Low Air Loss Mattress System

| Model # |  Dimensions |  Weight Capacity |  Flow Rate (Pump) |  Therapy Type |  Alarms |  Limited Warranty |
|---------|--|---|--|--|---|--|
| MA1000 | 80" L 36" W 10" H | 600 lb. | 1275 LPM | Alternating Pressure Pulsation True Low Air Loss | Audio & Visual | 1 Year (Cover/ Mattress) 1 Years (Power Unit) |
| MA1000M | microAIR MA1000 Mattress only | | | | | |
| MA1000P | microAIR MA1000 Pump only | | | | | |

***HCPCS Code: E0277**

©2019

Invacare Corporation.
All rights reserved.
Trademarks are identified
by the symbols ™ and
®. All trademarks are
owned by or licensed
to Invacare Corporation
unless otherwise noted.
Specifications are subject
to change without
notification.
Form No. 19-030
Rev. 12/19
191216

*All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.



Yes, you can.®