



INVACARE® MICROAIR® MA800
THERAPEUTIC SUPPORT SURFACE



Yes, you can.®

microAIR®

Invacare® microAIR® MA800

Alternating Pressure Low Air Loss Mattress System

**Medium - High Risk
Stage I-IV**

3-in-1 Alternating Pressure combined with Pulsation and True Low Air is Loss ideal for situations where maceration and moisture management are an issue



Highly vapor permeable and oversized pliable quilted nylon top cover designed to provide low shear, friction and moisture protection



Cable management system keeps cords away from the ground to help prevent trip hazards



Nylon non-slip base cover designed to help reduce shear and friction during bed articulation





Features & Benefits

10 adjustable comfort settings depending on the patient's therapy requirement

60 second pulsation feature is designed to help improve capillary blood flow by continually readjusting to a level surface pressure

Auto Firm Mode designed to provide quick maximum inflation within seconds to aid with transfers and nursing procedures

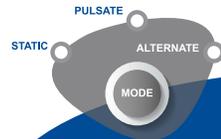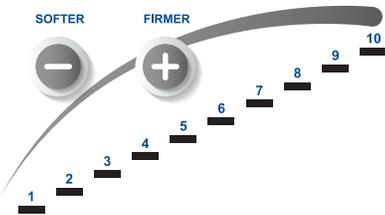

Advanced 3-in-1 alternating function inflates and deflates cells in a 3:1 cycle, providing active pressure relief by supporting 2/3 of the body at any given time

Customizable cycle times to meet patient's varying needs

True Low Air Loss Therapy provides extra comfort for moisture and temperature control




Foam base designed to help prevent bottoming out in the event of a power failure

Soft nylon urethane backed cells designed to provide easy cleaning and helps minimize the proliferation of bacterial growth

External exhaust keeps the mattress surface cooler and the pump quieter

HeelSense™ Technology designed to help provide further therapy and comfort by decreasing pressure in patient's vulnerable heel area

Request a Quote or Place an Order Today!

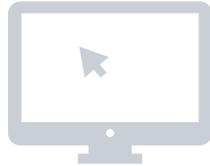

Contact your Inside Sales Representative

or



Contact Customer Service at 800.333.6900

Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2019

Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols ™ and ©. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.
Form No. 19-026
Rev. 12/19
191214

microAIR MA800 Alternating Pressure Low Air Loss Mattress System



Model #	Dimensions	Weight Capacity	Flow Rate (Pump)	Therapy Type	Alarms	Limited Warranty
MA800	80", 84" L 36", 42", 48" W 10" H	36" 600 lbs. 42", 48" 1000 lbs.	1275 LPM	Alternating Pressure Pulsation True Low Air Loss	Audio & Visual	1 Year (Cover/Mattress) 1 Year (Power Unit)
MA800	microAIR MA800 36" Alternating Pressure Low Air Loss Mattress System					
MA800B42	microAIR MA800 42" Alternating Pressure Low Air Loss Mattress System					
MA800B4284	microAIR MA800 42"W x 84"L Alternating Pressure Low Air Loss Mattress System					
MA800B48	microAIR MA800 48" Alternating Pressure Low Air Loss Mattress System					
MA800M	microAIR MA800 36" Mattress only					
MA800M42	microAIR MA800 42" Mattress only					
MA800M4284	microAIR MA800 42"W x 84" L Mattress only					
MA800M48	microAIR MA800 48" Mattress only					
MA800P	microAIR MA800 Pump only					

***HCPCS Code: E0277**

*All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.



Yes, you can.®